


## Analysis of Parabenes

Data No. LA104-0000

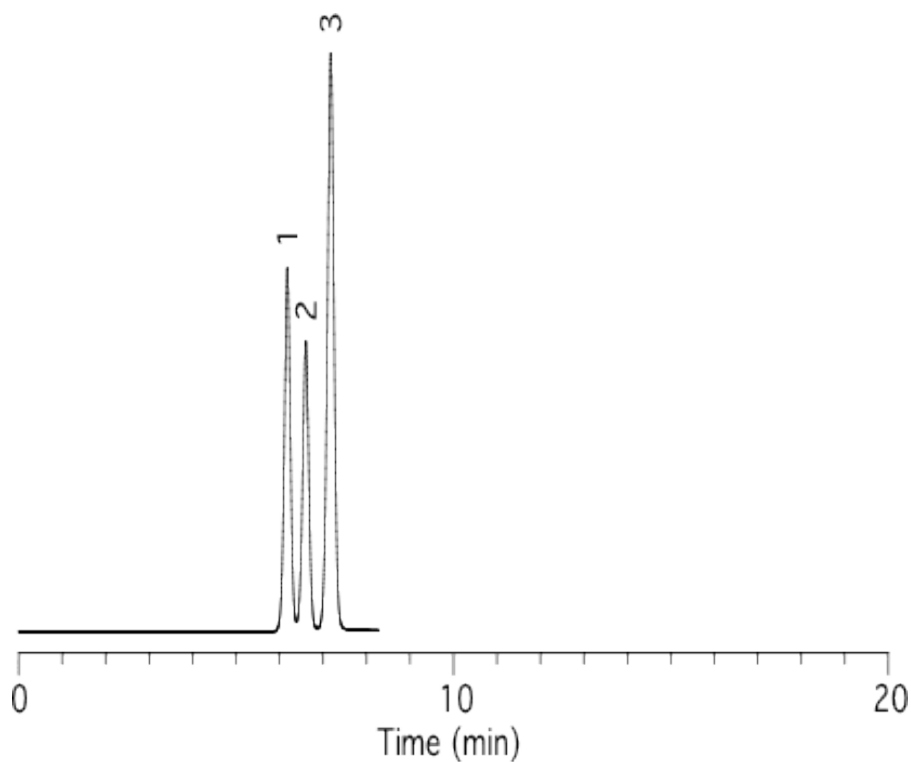

Column : Inertsil SIL-100A  
5  $\mu$ m  
150 x 4.6 mm I.D.

Solutes : 1. n-Butyl paraben  
2. n-Propyl paraben  
3. Ethyl paraben

Column Cat. No. : 5020-01711

Eluent : A) n-Hexane  
B) Ethanol  
A/B = 95/5, w/w

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :