


Free fatty acids

Data No. LA107-0000

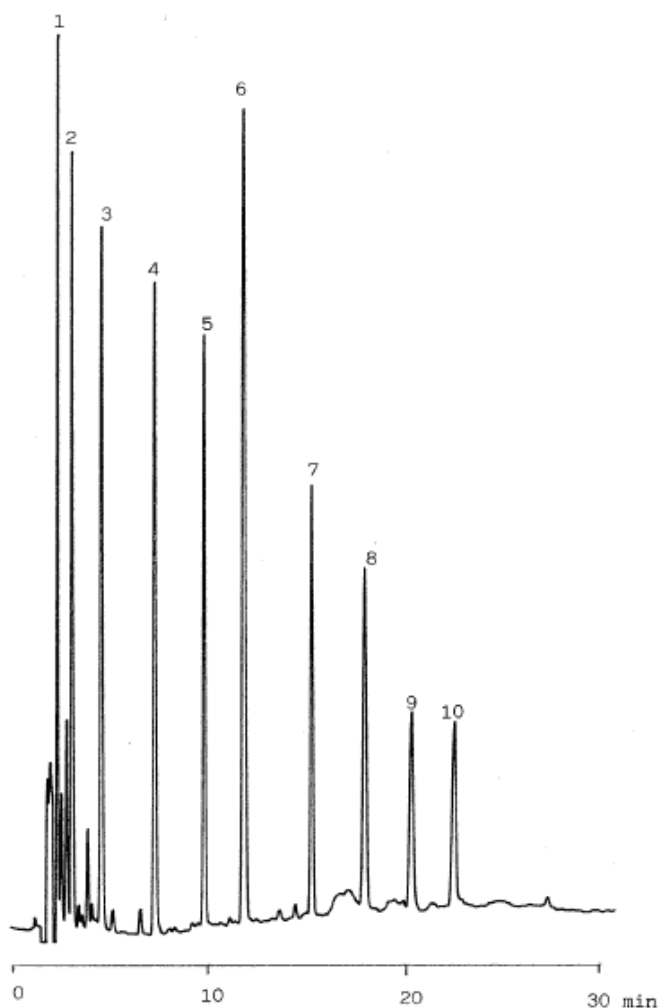

Solutes : 1. Propionic acid
2. Butyric acid
3. Valeric acid
4. Caproic acid
5. Heptylic acid
6. Caprylic acid
7. Capric acid
8. Lauric acid
9. Myristic acid
10. Palmitic acid

Column : Inertsil Ph[®]5 μ m.(150 \times 4.6mm I.D.)

Column Cat. No. : 5020-01301

Eluent : A:CH₃CN/0.1% H₃PO₄=25/75

B:CH₃CN/0.1% H₃PO₄=75/25

A/B? 100/0/20*min/0/100

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 210nm