


## Anion surfactants

Data No. LA111-0000

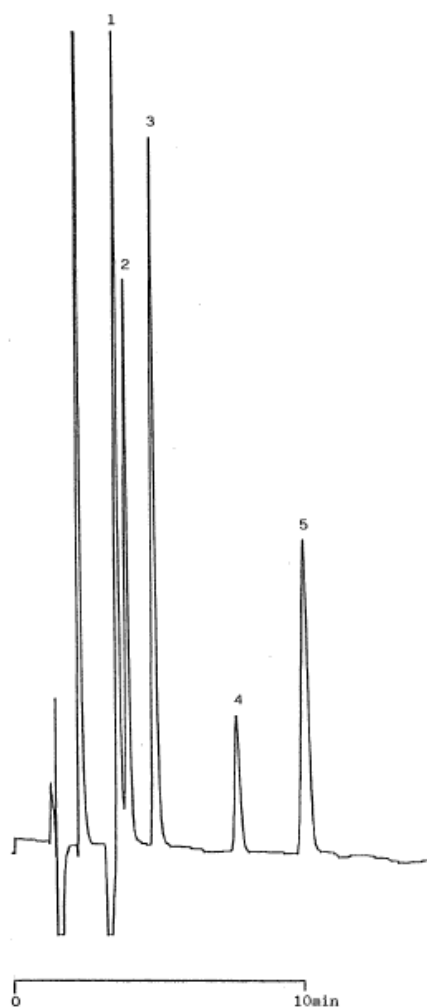

Column :Inertsil ODS-2  
5 $\mu$ m  
150  $\times$  4.6 mmI.D.

Column Cat. No. :5020-01101

Eluent :A:CH<sub>3</sub>CN/0.1M NaC<sub>10</sub>H<sub>7</sub>O<sub>4</sub>=80/20

Flow rate :1.0 mL/min.

Col. Temp. :40°C

Detection :RI

Injection Vol. :

Sample :

Solutes : 1. n-Undecane sulfonate Na  
2. n-Dodecane sulfonate Na  
3. n-Tridecane sulfonate Na  
4. n-Tetradecane sulfonate Na  
5. n-Hexadecane sulfonate Na