


Anion surfactants

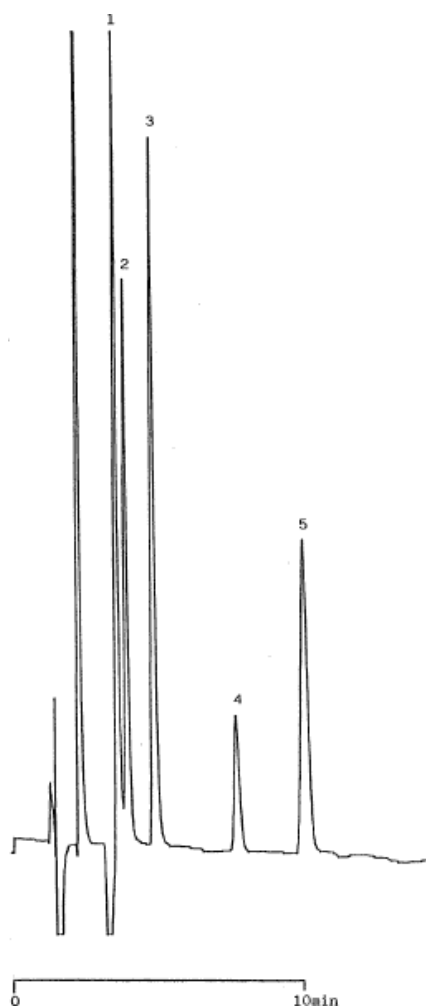

Column : Inertsil ODS-2
 5 μ m
 150 \times 4.6 mm I.D.

Column Cat. No. : 5020-01101

Eluent : A: CH₃CN/0.1M NaClO₄=80/20

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : RI

Injection Vol. :

Sample :

Solutes : 1. n-Undecane sulfonate Na
 2. n-Dodecane sulfonate Na
 3. n-Tridecane sulfonate Na
 4. n-Tetradecane sulfonate Na
 5. n-Hexadecane sulfonate Na