


InertSearch™ for LC

Inertsil® Applications

Organic acids

Data No.LA120-0000

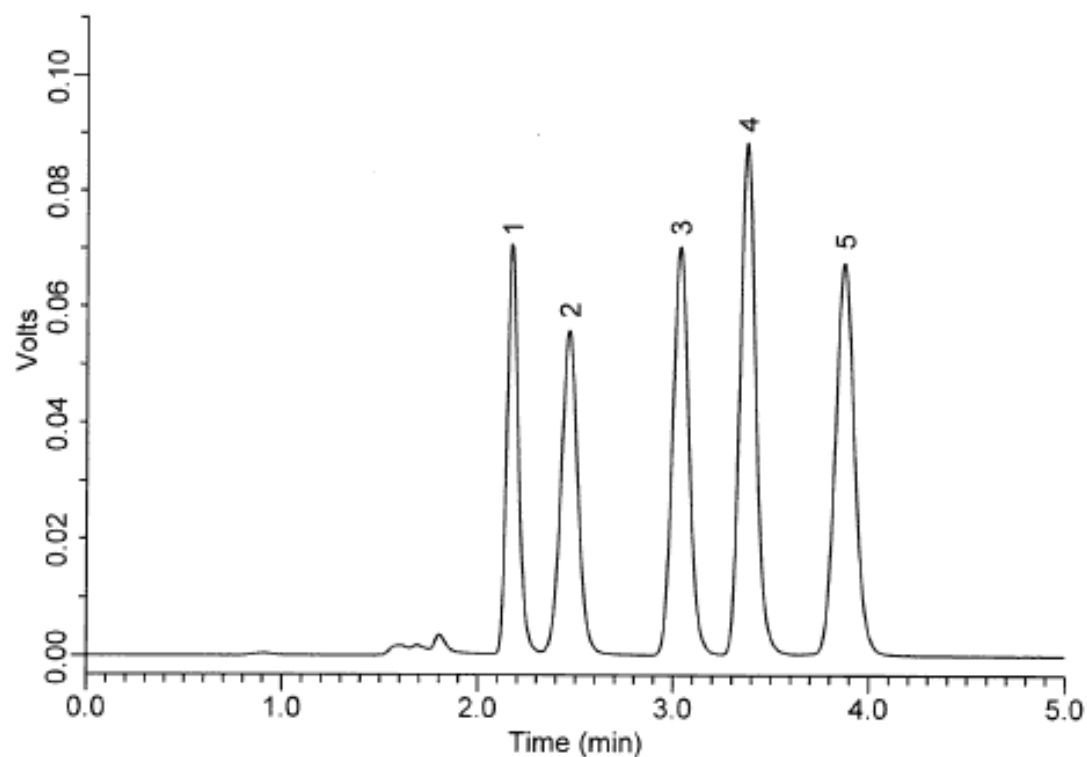

Column	: Inertsil C8-3 150 × 4.6mm I.D.	Solutes	: 1. Formic acid(0.5mg/mL) 2. alpha-Ketoglutaric acid(0.05mg/mL) 3. Citric acid(0.5mg/mL) 4. Maleic acid(0.005mg/mL) 5. Fumaric acid(0.005mg/mL)
Column Cat. No. :	5020-01900		
Eluent	: (NH ₄)H ₂ PO ₄ (pH 2.8; H ₃ PO ₄)		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 210 nm		
Injection Vol.	:		
Sample	: Organic acids		