


Anion Surfactants

Data No. LA123-0000

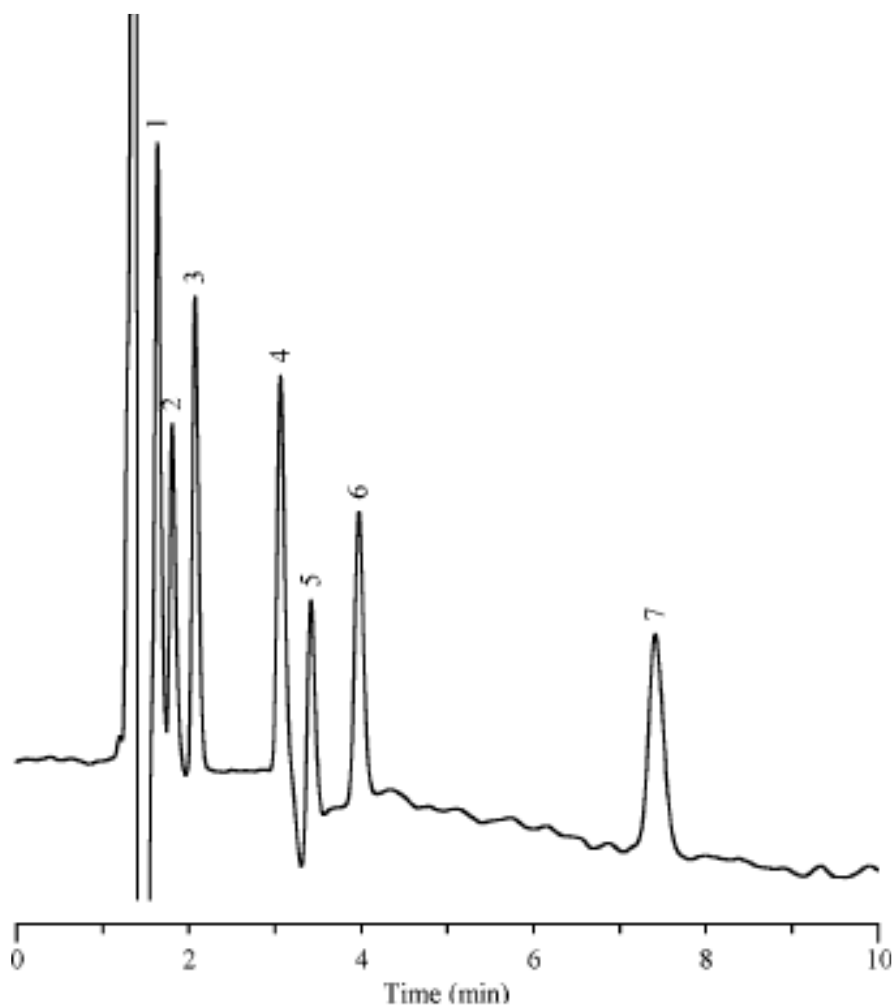

Column : Inertsil ODS-3V
5µm
150×4.6mmI.D.

Column Cat. No. : 5020-01801

Eluent : 0.1M NaClO₄ / acetonitrile (40/60)

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : RI

Injection Vol. :

Sample : Anion Surfactants (80mg/mL)

Solutes :1. Nonyl sulfonate Na

2. Decyl sulfate Na

3. Undecyl sulfonate Na

4. Tridecyl sulfonate Na

5. Dodecyl sulfate Sodium

7. Tetradecyl sulfonate Na

8. Hexadecyl sulfonate Na