


InertSearch™ for LC

Inertsil® Applications

Tartaric acid

Data No. LA124-0000

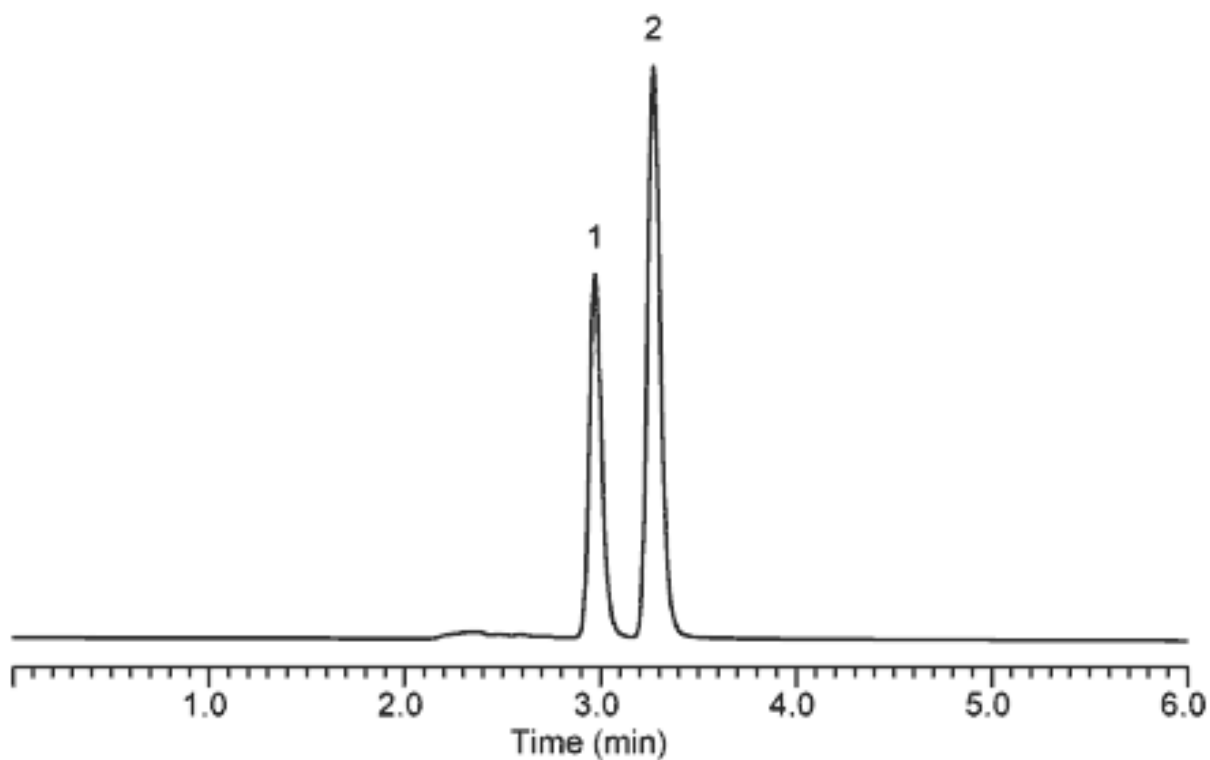

Column : Inertsil ODS-3V
5µm
250×4.6mmI.D.

Solutes : 1. meso-Tartaric acid 2.970 237144 8802.83
2. DL-Tartaric acid 3.270 396050 10167.9

Column Cat. No. : 5020-01802

Eluent : 50mM Phosphate buffer (pH2.5)

Flow rate : 1mL/min

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 5µL

Sample : meso-Tartaric acid 1-hydrate (0.5 mg/mL)
DL-Tartaric acid (0.5 mg/mL)