


InertSearch™ for LC

Inertsil® Applications

Phthalic acid

Data No. LA125-0000

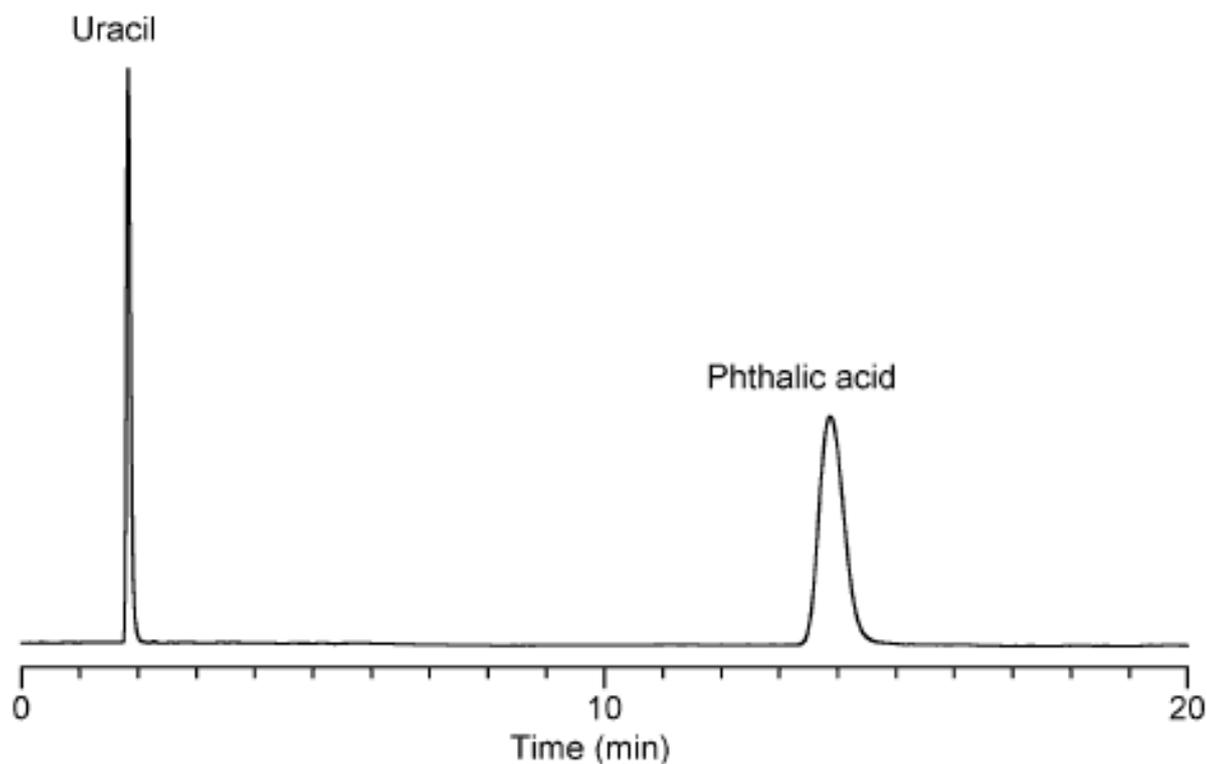

Column	: Inertsil ODS-3V 5 μ m 150 \times 4.6mmI.D.	Solutes	: 1. Uracil 1.830 143196 2. Phthalic acid 13.870 334340
Column Cat. No.	: 5020-01801		
Eluent	: Acetonitrile / 0.1% H ₃ PO ₄ = 10 / 90		
Flow rate	: 1.0mL/min		
Col. Temp.	: 40°C		
Detection	: UV 210 nm		
Injection Vol.	: 5 μ L		
Sample	: Uracil Phthalic acid		