


Organic acids

Data No. LA127-0000

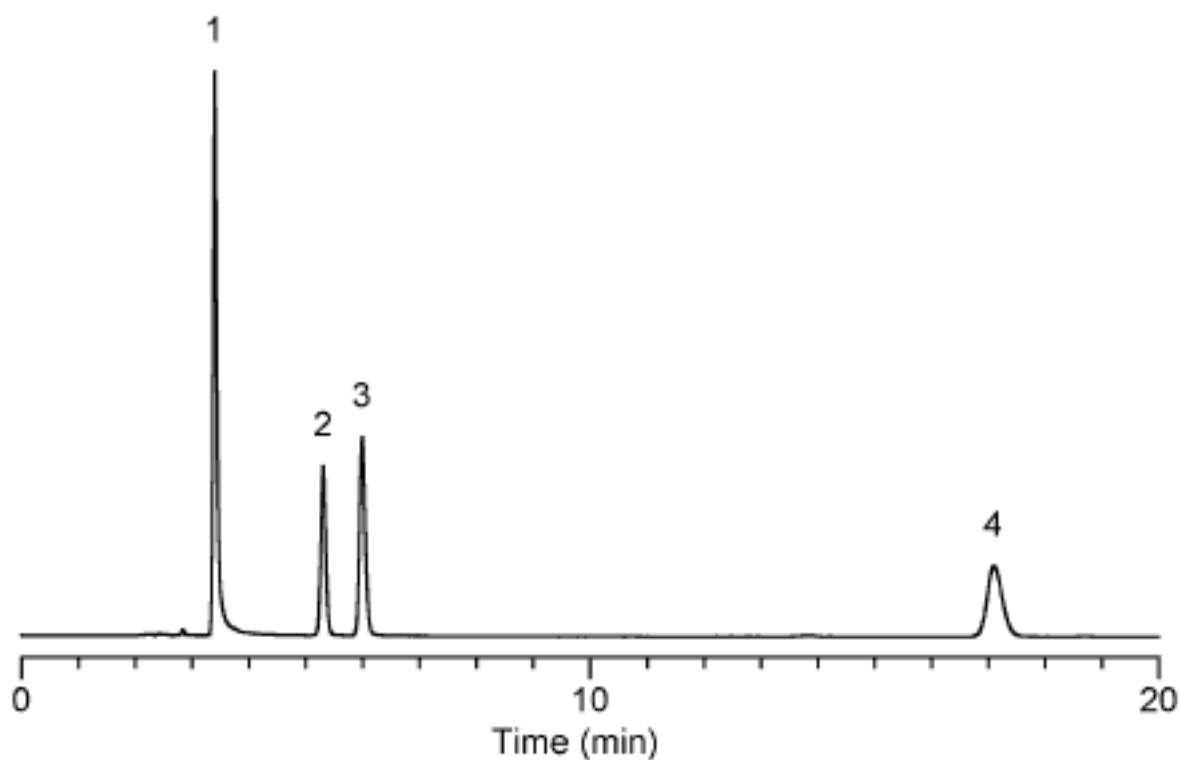

Column : Inertsil ODS-3V
5 μ m
250 \times 4.6mmI.D.

Column Cat. No. : 5020-01802

Eluent : 20mMKH₂PO₄ in H₂O (pH 3.0, H₃PO₄)

Flow rate : 1.0mL/min

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 20 μ L

Sample : Organic acids

Solutes : 1. Formic acid (0.1 %)
2. Lactic acid (0.1 %)
3. Acetic acid (0.1 %)
4. Propionic acid (0.1 %)