


InertSearch™ for LC

Inertsil® Applications

Diphenhydramine in Medicine

Data No. LA128-0000

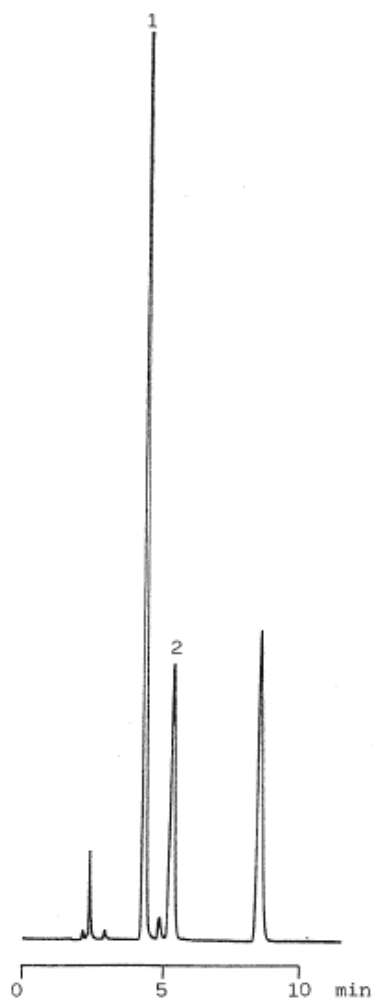

Column : Inertsil C8
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. Methylsalicylate
2. Diphenhydramine

Column Cat. No. : 5020-01201

Eluent : CH₃CN/10mM IPCC-12A=55/45
IPCC-12A: Sodium 1-Dodecylsulfate

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :