


InertSearch™ for LC

Inertsil® Applications

Hinokitiol in Cosmetic

Data No. LA130-0000

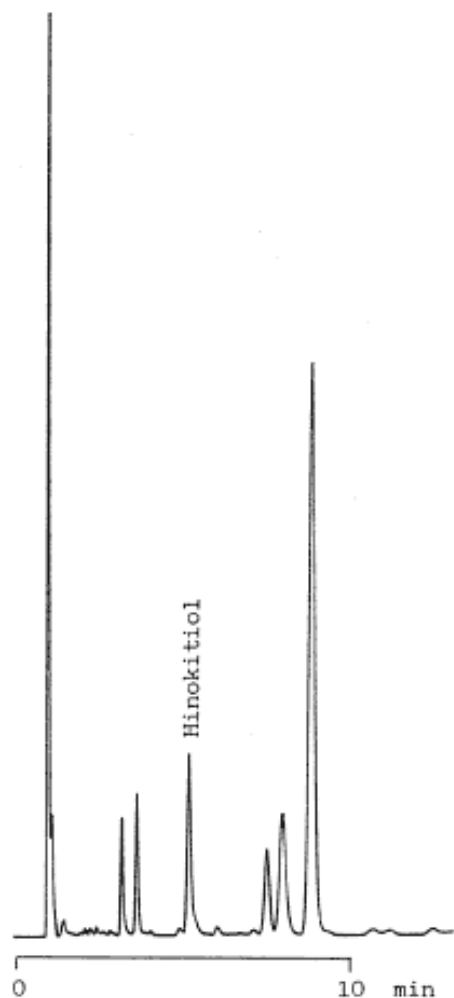

Column : Inertsil C4

5 μ m

150 \times 4.6mm I.D.

Column Cat. No. : 5020-01401

Eluent : CH₃CN/10mM KH₂PO₄+0.1%EDTA-2Na
(pH3.5 by H₃PO₄)=35/65

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample : Cosmetic

Solutes : Hinokitiol