


Herb medicine

Data No. LA133-0000

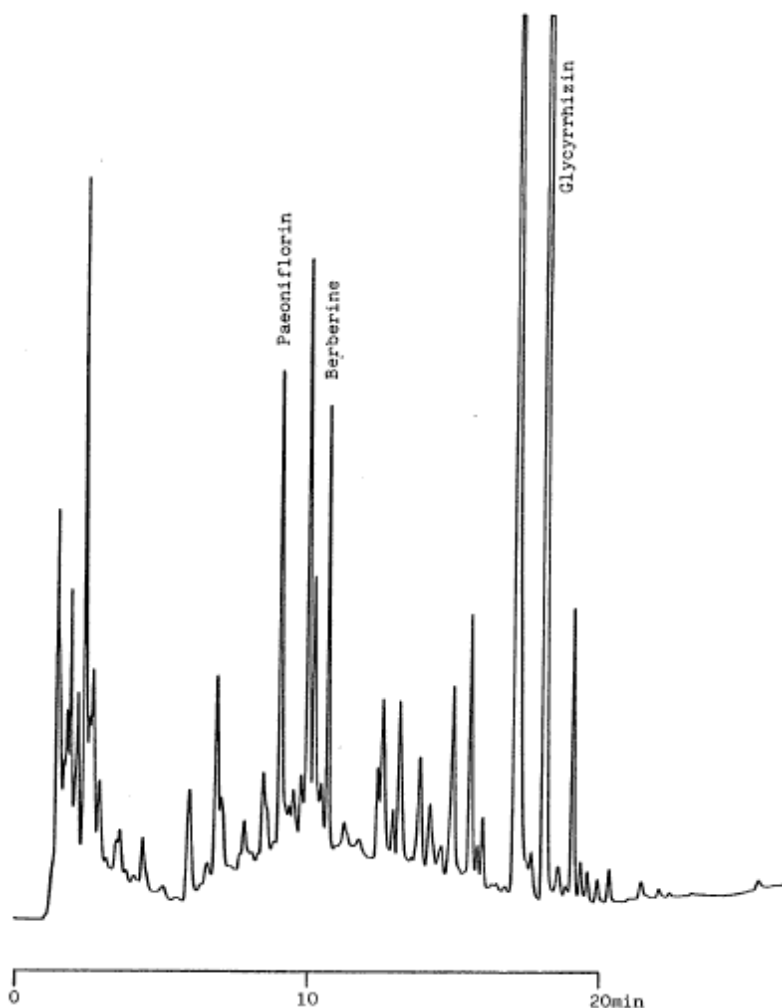

Column : Inertsil ODS-2
5 μ m
150 \times 4.6 mm I.D.
Column Cat. No. : 5020-01101
Eluent : A: CH₃CN
B: 0.05% H₃PO₄
A/B 10/90 \rightarrow (25min) \rightarrow 90/10
Flow rate : 1.0 mL/min.
Col. Temp. : 40°C
Detection : UV 254 nm
Injection Vol. :
Sample :
Solutes : Paoniflorin
Berberin
Glycyrrhizin