


Scutellariae radix

Data No. LA135-0000

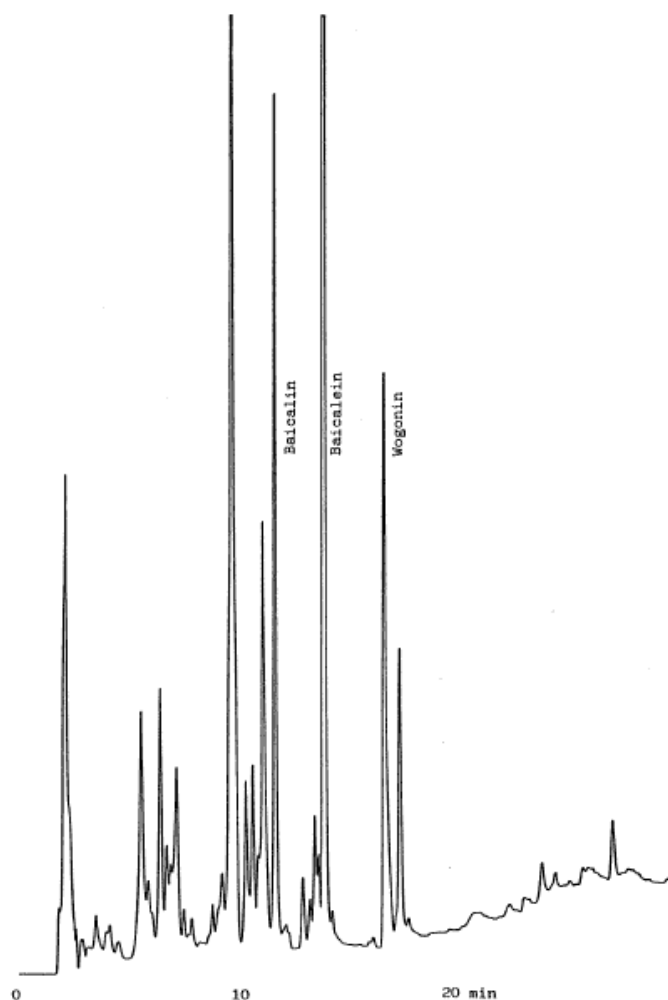

Column : Inertsil ODS-2
5 μ m
150 \times 4.6 mmI.D.

Column Cat. No. : 5020-01102

Eluent : A: CH₃CN
B: 0.05% CF₃COOH
A/B 20/80 \rightarrow (20min) \rightarrow 80/20

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Baicalin
2. Baicalein
3. Wogonin