


InertSearch™ for LC

Inertsil® Applications

Insulin standards E10

Data No. LA137-0000

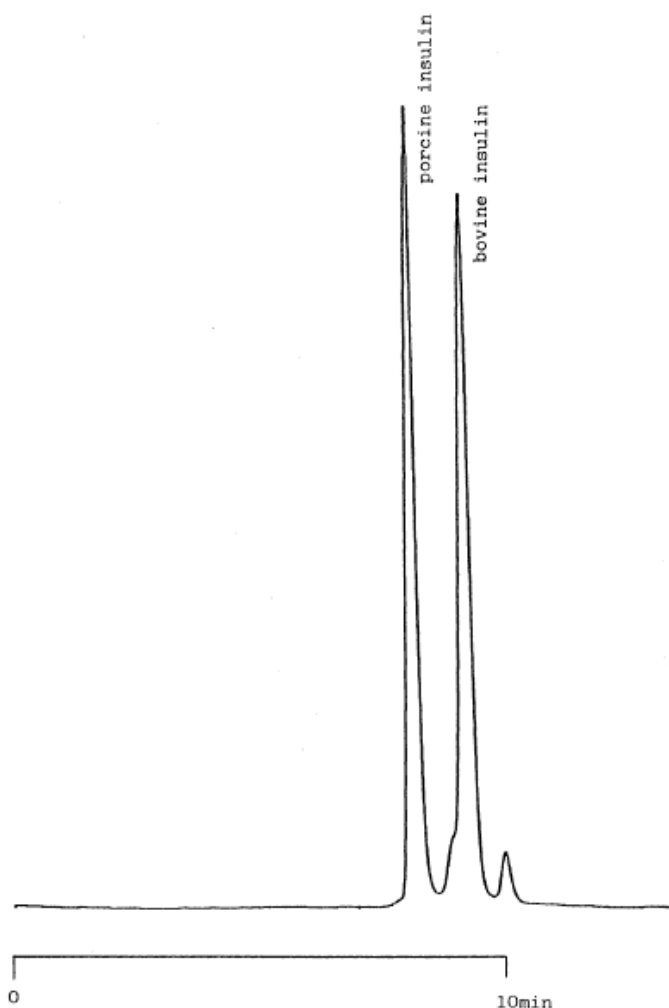

Column : Inertsil ODS-2
5 μ m
250 \times 4.6 mmI.D.
Column Cat. No. : 5020-01102
Eluent : A:0.1% CF₃COOH in CH₃CN
B:0.1% CF₃COOH
A/B 20/80?(5min)?30/70
→(40min)→50/50
Flow rate : 1.0 mL/min.
Col. Temp. : RT
Detection : UV 280 nm
Injection Vol. :
Sample :
Solutes :porcine insulin
Bovine insulin