


InertSearch for LC Data No. LA138-0000

Components of Cold-drugs

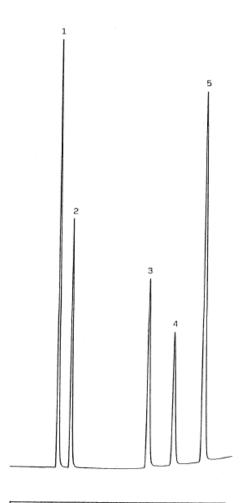

Column : Inertsil ODS-2

5µm

 150×4.6 mmI.D.

Column Cat. No. :5020-01101

Eluent : A:CH₃OH/0.1% H₃PO₄=30/70

B:CH3OH/0.1% H3PO4=80/20

10min

A/B 100/0→(15min)→0/100

Flow rate : 1.0 mL/min.

Col. Temp. :40°C

Detection : UV 220 nm

Injection Vol. : Sample : Solutes

: 1. Acetaminophen

2. Caffeine

3. Ethenzamide

4. Bucetin

5. Acetyl salicylic acid