


Purification of Glycyrrhizin

Data No. LA139-0000

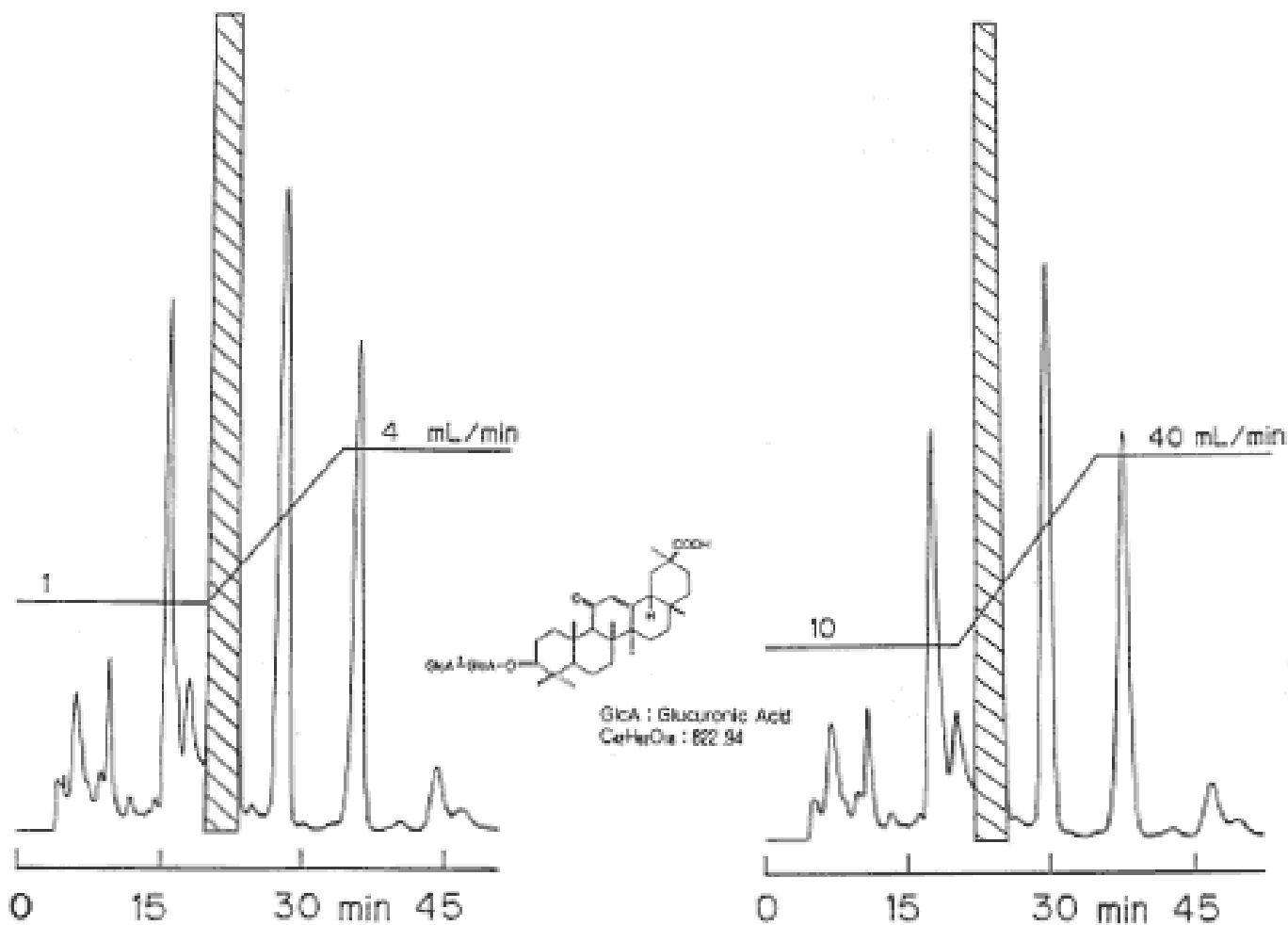

Column : Inertsil PREP-ODS
10µm

Column Cat. No. : 5020-35121

Eluent : CH₃CN/10% CH₃COOH=35/65

Flow rate : Flow gradient

Col. Temp. :

Detection : UV 254 nm

Injection Vol. :

Sample : crude Glycyrrhizin

250 × 6.0 mmI.D. - Column size - 250 × 20 mmI.D.

1.0→4.0mL/min - Flow rate - 10→40mL/min

300µL - Sample size - 3000µL

Solutes : 1. Glycyrrhizin