


Purification of Zingiberis rhizoma

Data No. LA140-0000

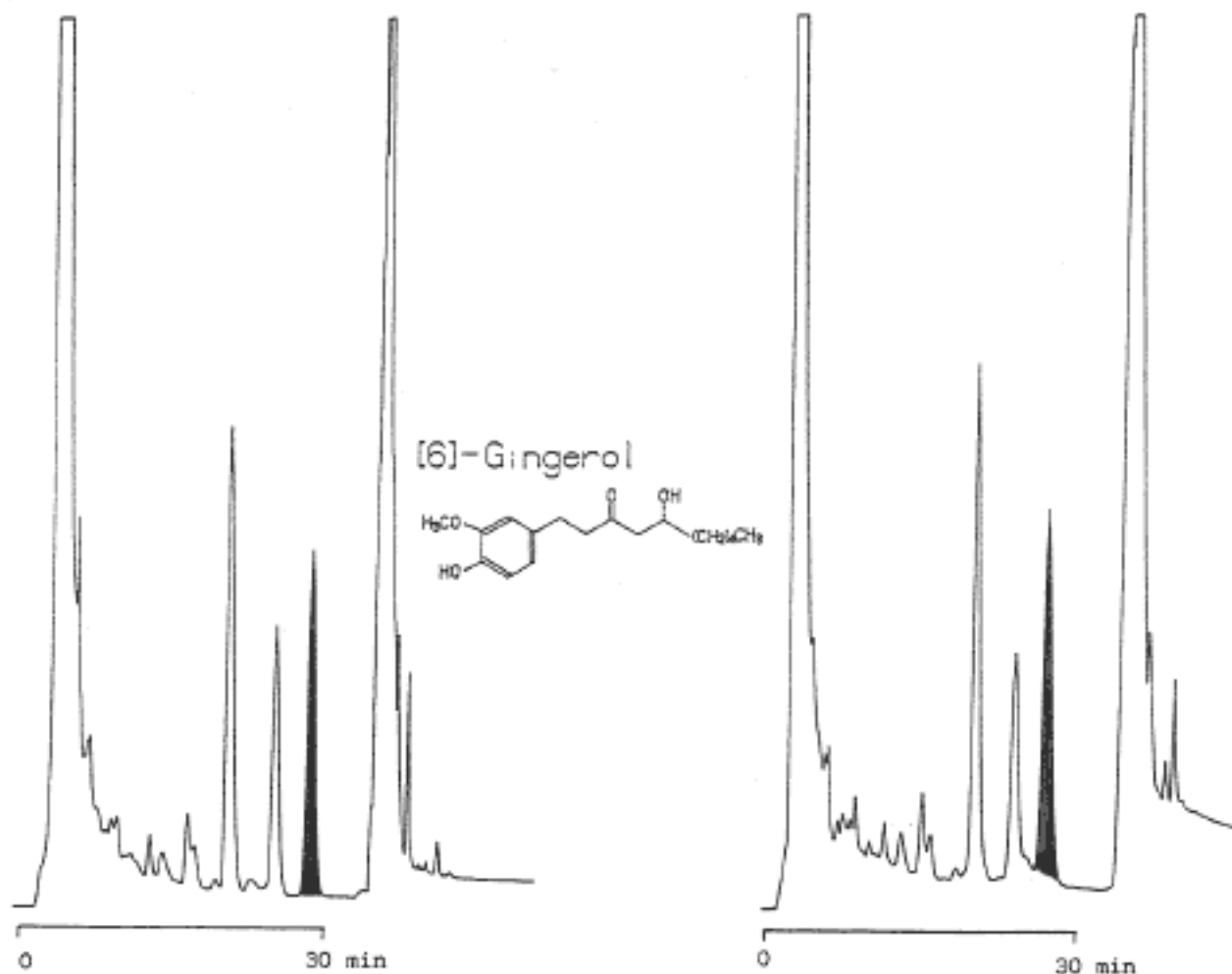

Column : Inertsil PREP-ODS
10 μ m

Column Cat. No. : 5020-35151

Eluent : A:CH₃OH

B:H₂O

A/B 60/40→(35min)→100/0

Flow rate :

Col. Temp. :

Detection : UV 254 nm

Injection Vol. :

Sample : Zingiberis rhizoma

250 × 6.0 mmI.D. - Column size - 250 × 50 mmI.D.

1mL/min - Flow rate - 70mL/min

0.5mL - Sample size - 35mL

Solutes : 1. Zingiberis rhizoma