


Recycle preparative purification of Steroids

Data No. LA141-0000

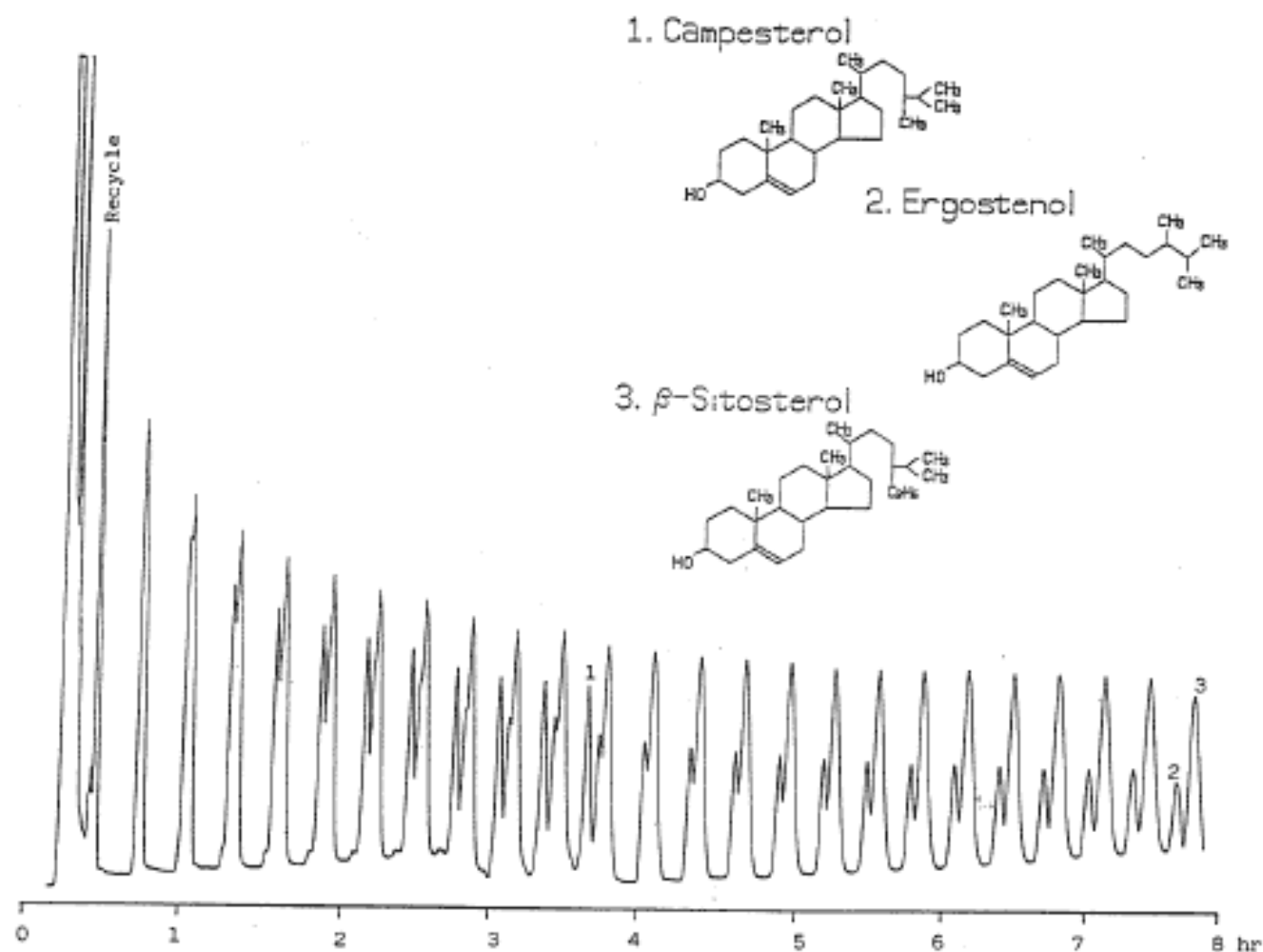

Column : Inertsil PREP-ODS
10 μ m
250 \times 20mm I.D.

Solutes : 1. Campesterol
2. Ergosterol
3. β -Sitosterol

Column Cat. No. : 5020-31123

Eluent : C₂H₅OH/H₂O=93/7

Flow rate : 12 mL/min.

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 2mL

Sample : Crude β -sitosterol(20mg/mL)