


## Purification of Puerarie radix

Data No. AL142-0000

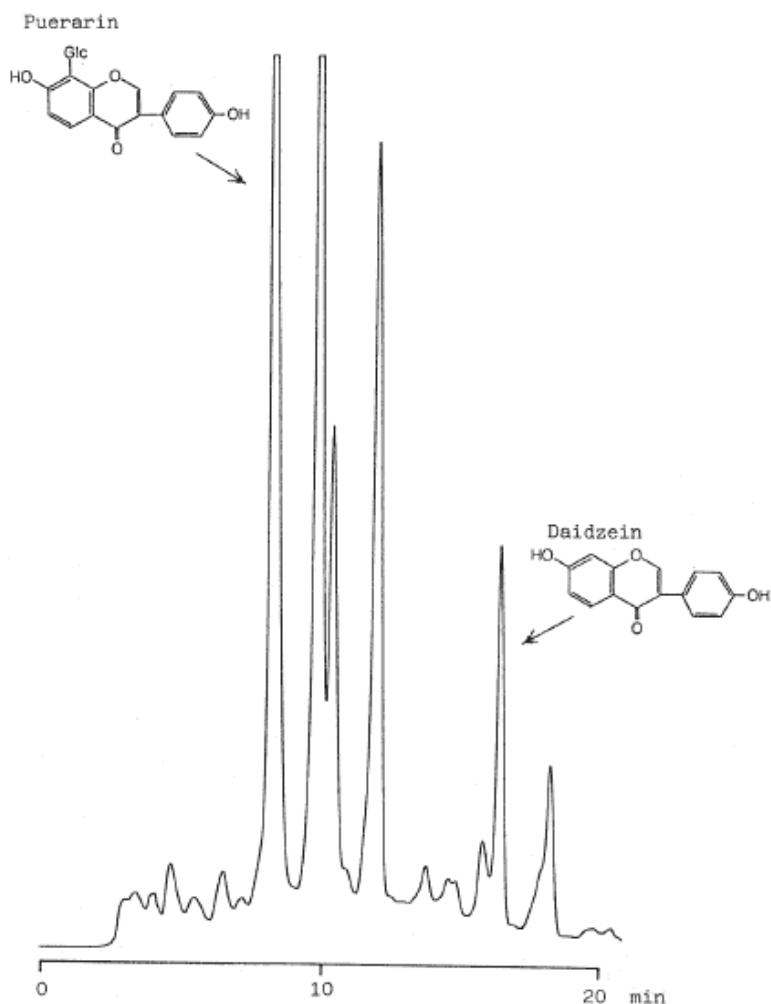

Column : Inertsil PREP-ODS  
10µm  
250 × 20mmI.D.

Solutes : 1. Puerarin  
2. Daidzein

Column Cat. No. : 5020-31123

Eluent : A:CH<sub>3</sub>OH  
B:H<sub>2</sub>O  
A/B 30/70→(20min)→90/10

Flow rate : 15 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. : 10mL

Sample :