


## Sulphur drugs

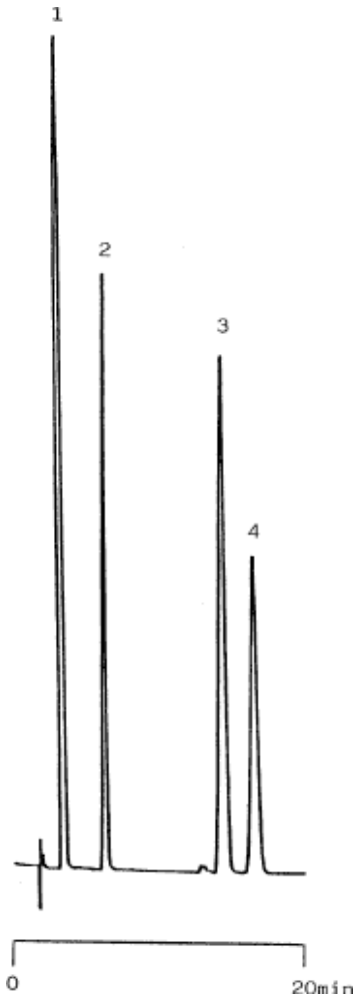

Column : Inertsil ODS-2  
5 $\mu$ m  
150  $\times$  4.6mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH<sub>3</sub>OH/5mM IPCC-07=28/72  
IPCC-07:Sodium 1-Heptanesulfonate

Flow rate : 1.2 mL/min.

Col. Temp. : RT

Detection : UV 254 nm

Injection Vol. :

Sample :

Solutes : 1. Sulfathiazole  
2. Sulfamethizole  
3. Sulfadimethoxine  
4. Sulfaquinoxaline