


Purification of *Coptidis rhizoma*

Data No. LA145-0000

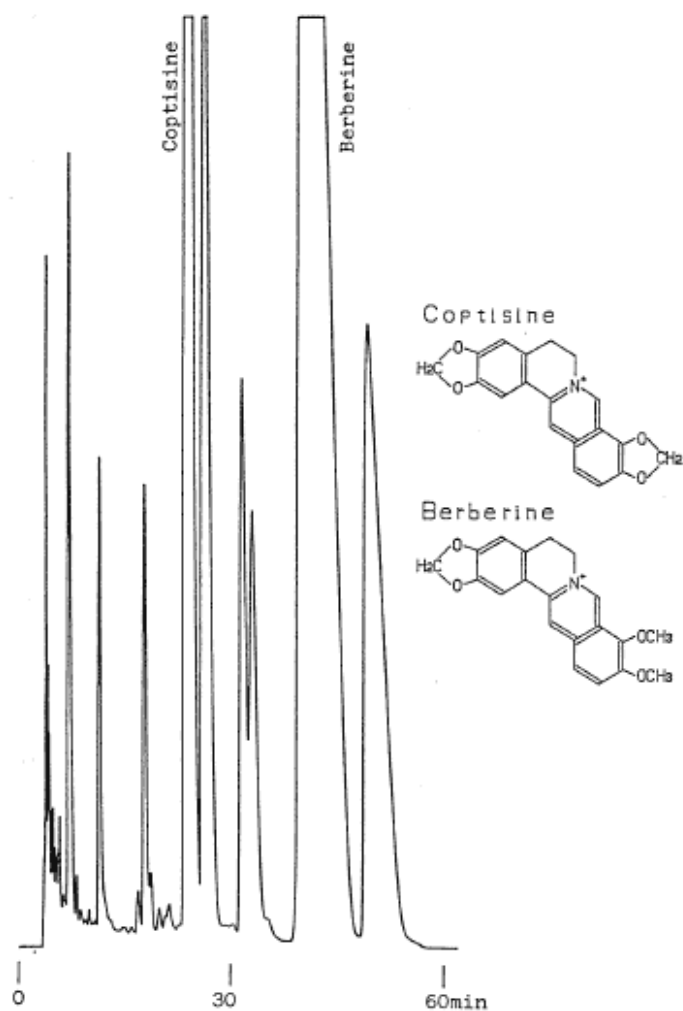

Column : Inertsil PREP-ODS
10µm
250 × 30mmI.D.

Solutes : 1. Coptisine
2. Berberine

Column Cat. No. : 5020-31133

Eluent : CH₃OH/0.1% CF₃COOH=30/70

Flow rate : 30 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. : 2.5mL

Sample : *Coptidis Rhizoma*(50mg/mL)