


Purification of Rhei rhizoma

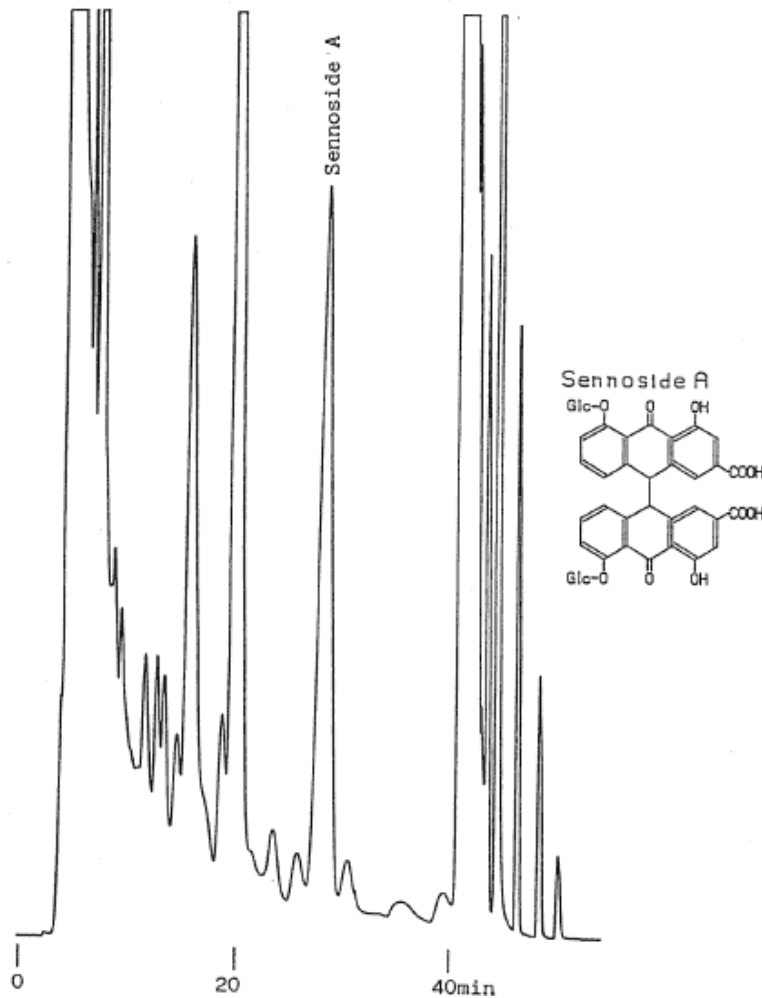

Column	: Inertsil PREP-ODS 10 μ m 250 \times 20mmI.D.	Solutes	: 1. SennosideA
Column Cat. No.	: 5020-31123		
Eluent	: CH ₃ OH/0.1% CF ₃ COOH=38/62		
Flow rate	: 10 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 270 nm		
Injection Vol.	: 1. 2mL		
Sample	: Rhei Rhizoma(50mg/mL)		