


InertSearch™ for LC

Inertsil® Applications

Ipecacuanhae radix

Data No. LA147-0000

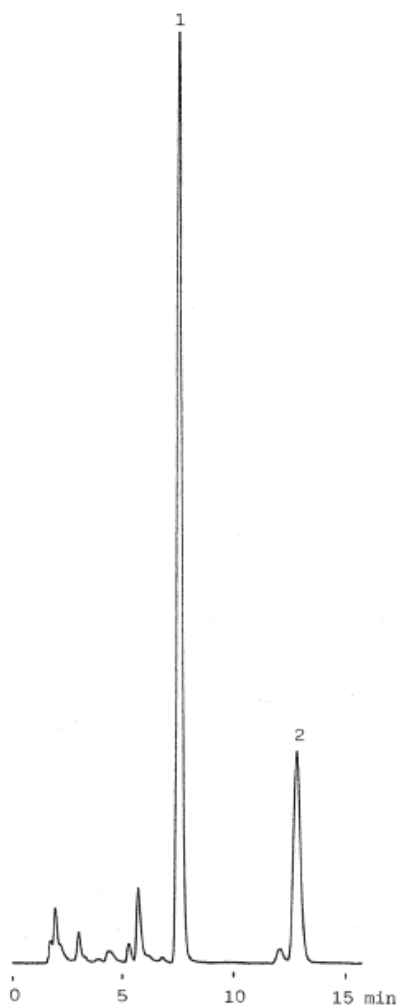

Column : Inertsil ODS-2
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃CN/5mM IPCC-07+0.1% H₃PO₄=25/75
IPCC-07:Sodium 1-Heptansulfonate

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : FL
Ex 283nm
Em 318nm

Injection Vol. :

Sample :

Solutes : 1. Cephaeline
2. Emetine