


InertSearch™ for LC

Inertsil® Applications

Butyrophenone derivatives

Data No. LA160-0000

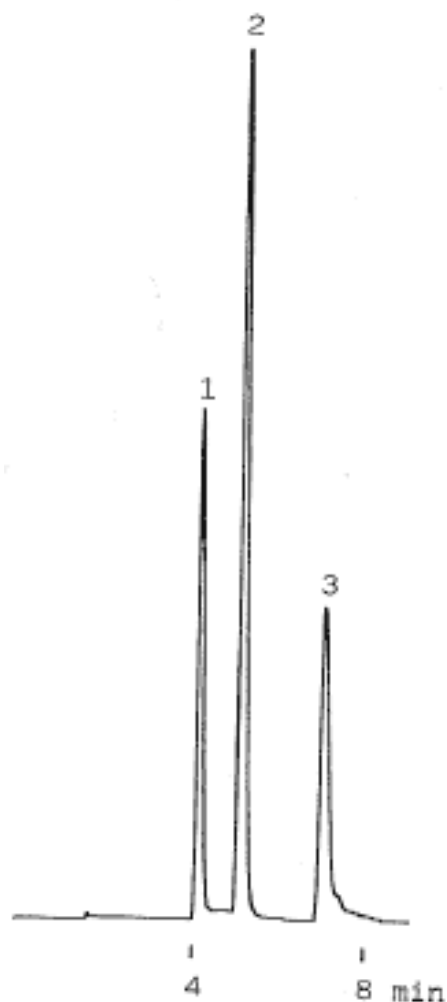

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. Droperidol
2. Spiroperidol
3. Haloperidol

Column Cat. No. : 5020-13001

Eluent : CH₃CN/50mM KH₂PO₄=35/65

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :