


Analysis of Tocopherol

Data No. LA163-0000

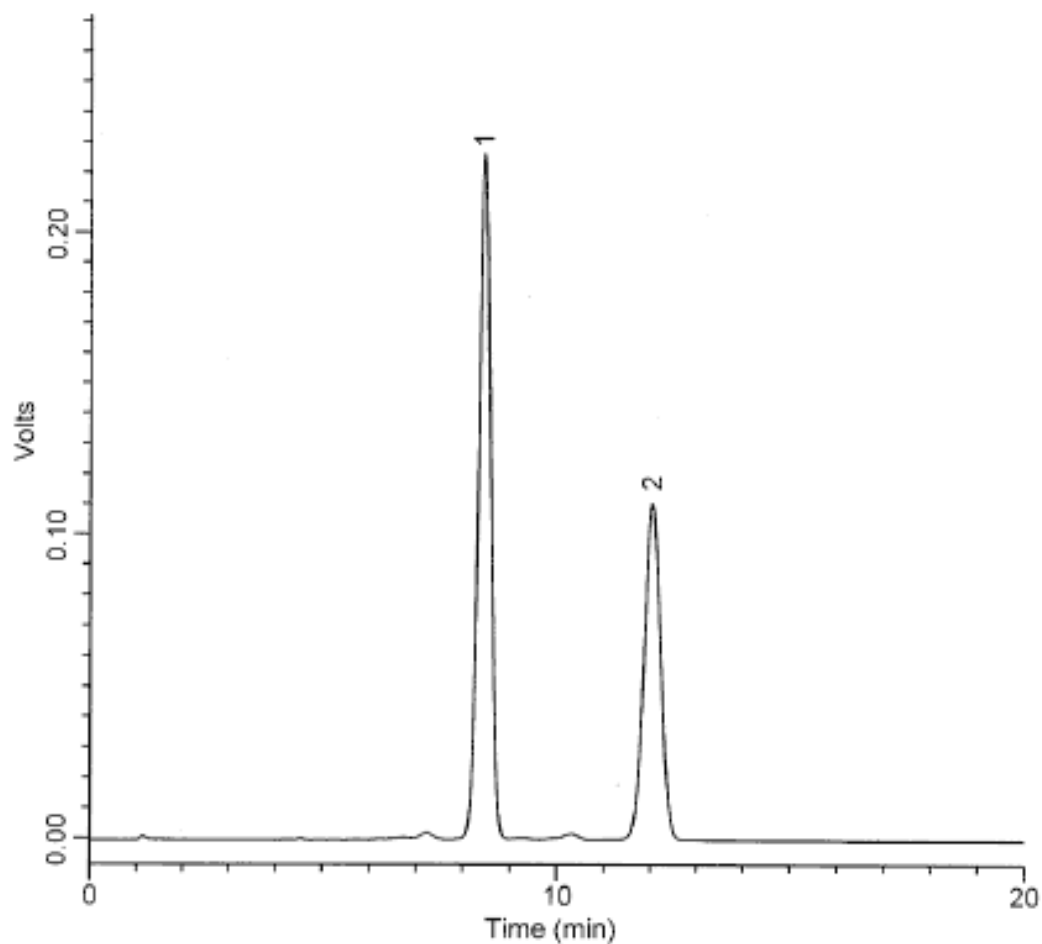

Rs=6.3

Column : Inertsil ODS-3V
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. DL- α -Tocopherol(1.0mg/mL)
2. DL- α -Tocopherol acetate(1.0mg/mL)

Column Cat. No. :

Eluent : CH₃OH/H₂O=98/2

Flow rate : 1.6 mL/min.

Col. Temp. : 30°C

Detection : UV 284 nm

Injection Vol. : 20 μ L

Sample : 1. DL- α -Tocopherol(1.0mg/mL)
2. DL- α -Tocopherol acetate(1.0mg/mL)