


Analysis of Procaterol hydrochloride

Data No. LA165-0000

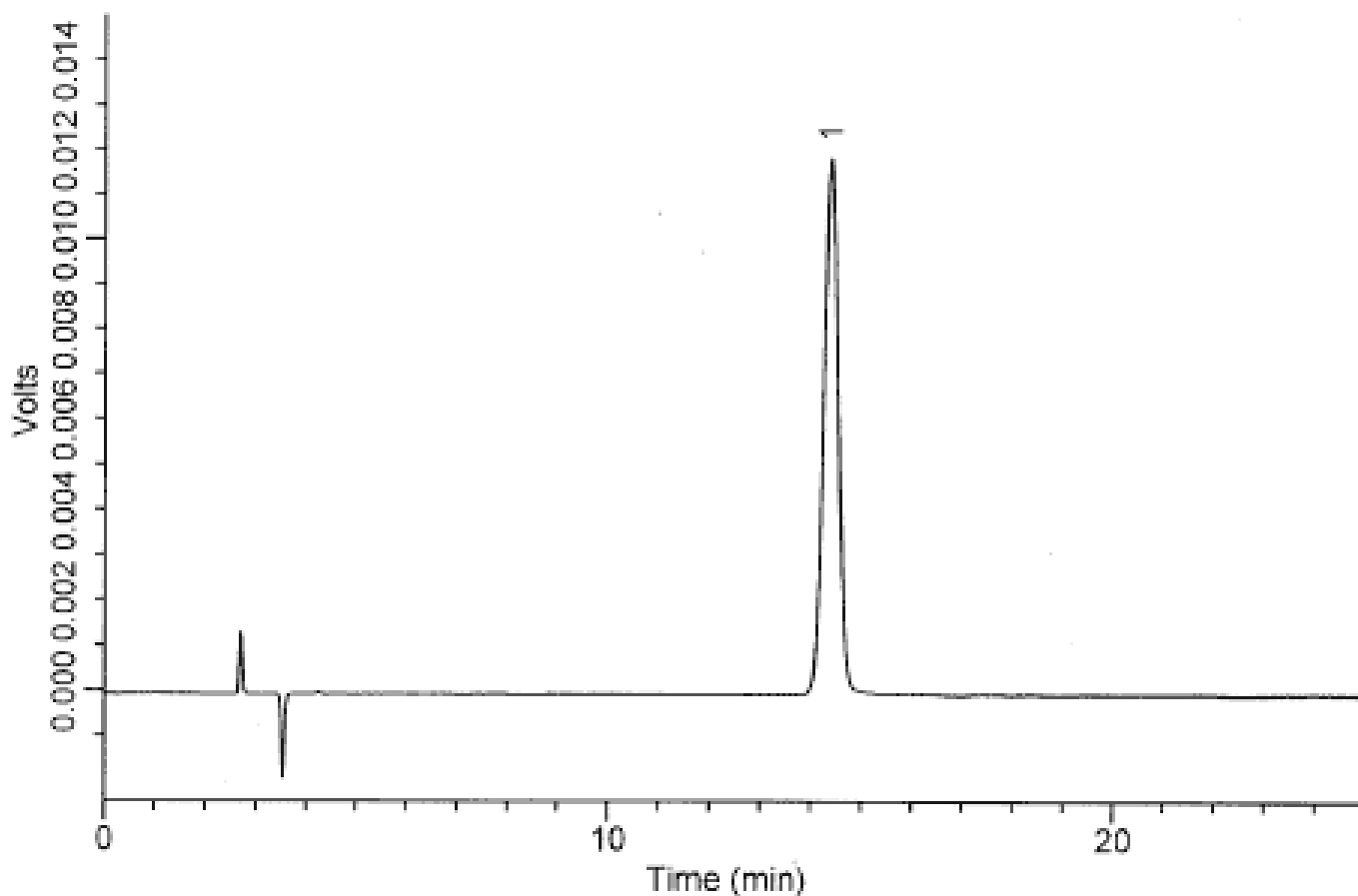

Column : Inertsil ODS-3V
250 × 4.6mm I.D. Solutes : 1. Procaterol hydrochloride(25µg/mL)

Column Cat. No. :

Eluent : CH₃OH/5mM IPC-5/CH₃COOH=23/76/1

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. : 2µL

Sample : 1. Procaterol hydrochloride(25µg/mL)