


Analysis of Fluocinolone acetonide

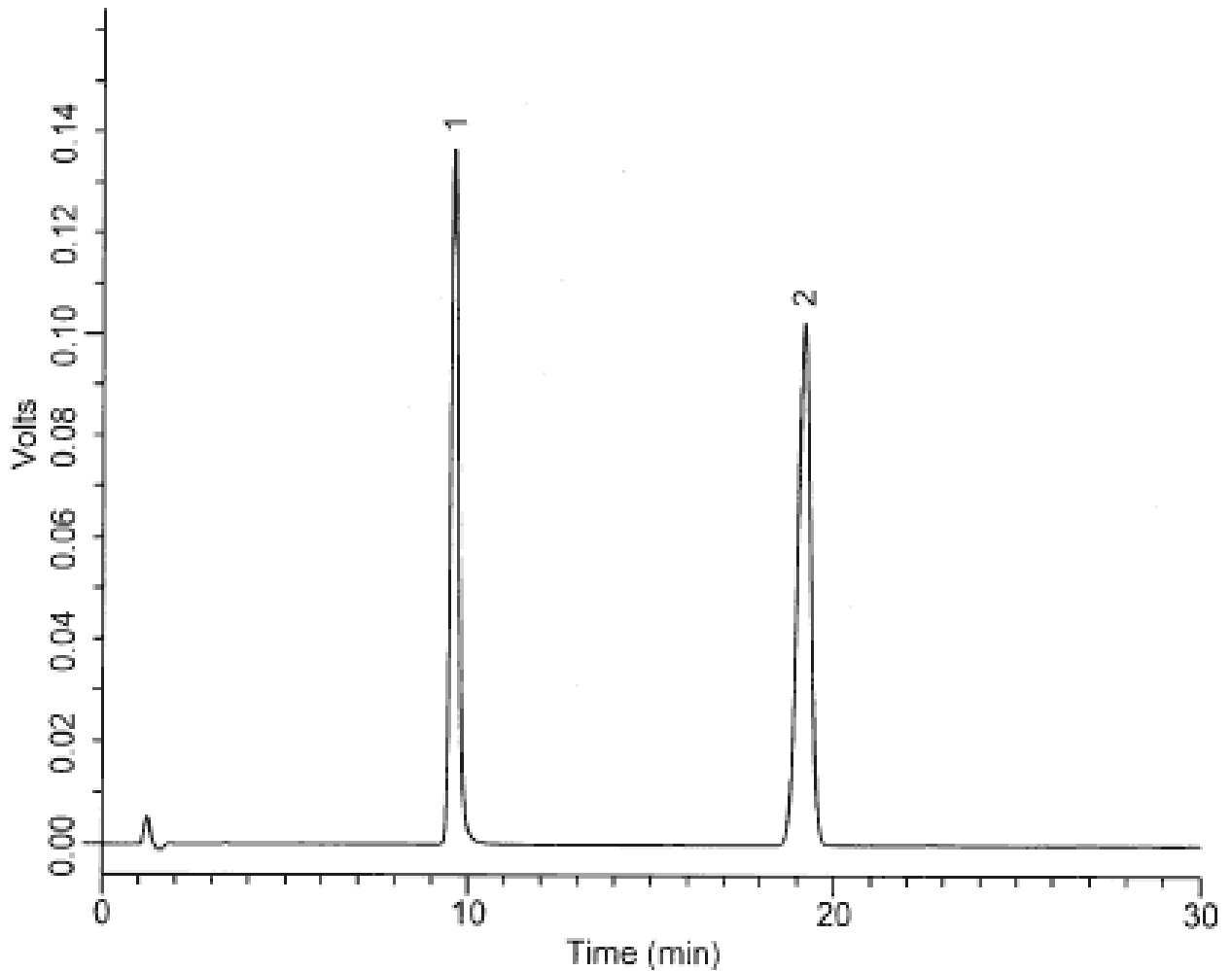

Column	:Inertsil ODS-3V 5µm 250 × 4.6mmI.D.	Solutes	: 1. Ethyl paraben(0.04mg/mL) 2. Fluocinolone acetonide(0.2mg/mL)
Column Cat. No.	:		
Eluent	:CH ₃ CN/H ₂ O=30/70		
Flow rate	:2.0 mL/min.		
Col. Temp.	:40°C		
Detection	:UV 254 nm		
Injection Vol.	:20µL		
Sample	:1. Ethyl paraben(0.04mg/mL) 2. Fluocinolone acetonide(0.2mg/mL)		