


Analysis of Fluocinolone acetonide

Data No. LA170-0000

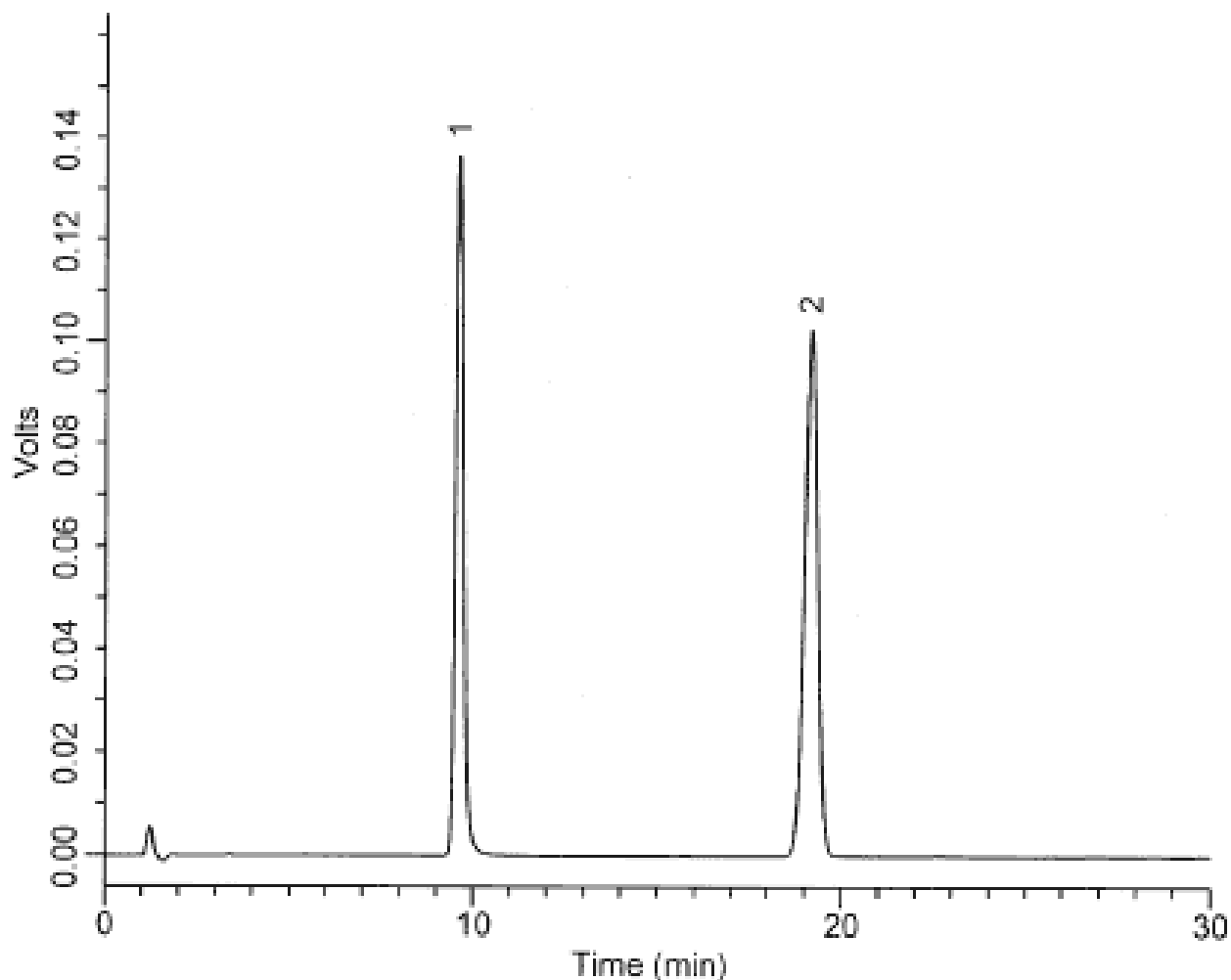

Column : Inertsil ODS-3V
5 μ m
250 \times 4.6mm I.D.

Solutes : 1. Ethyl paraben(0.04mg/mL)
2. Fluocinolone acetonide(0.2mg/mL)

Column Cat. No. :

Eluent : CH₃CN/H₂O=30/70

Flow rate : 2.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. : 20 μ L

Sample : 1. Ethyl paraben(0.04mg/mL)
2. Fluocinolone acetonide(0.2mg/mL)