


Analysis of Batamethasone valerate

Data No. LA174-0000

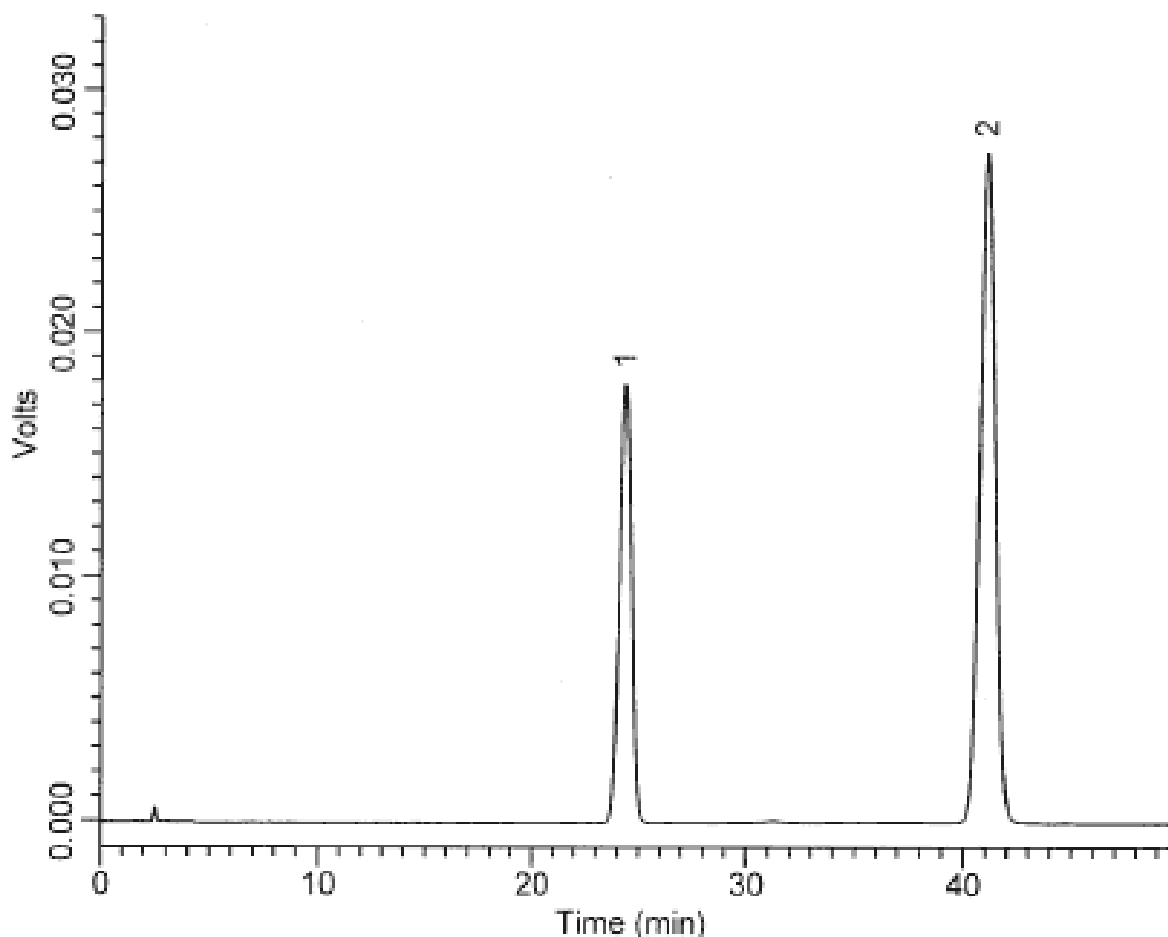

$R_s=14.2$

Column : Inertsil ODS-3V
5 μ m
250 \times 4.6mm I.D.

Solutes : 1. Betamethasone valerate (0.05 mg/mL)
2. Isoamyl benzoate (0.5 μ g/mL)

Column Cat. No. :

Eluent : CH₃OH/H₂O=70/30

Flow rate : 1.0 mL/min.

Col. Temp. : 25°C

Detection : UV 254 nm

Injection Vol. : 10 μ L

Sample : 1. β -methasone zalerate (0.05mg/mL)
2. Isoamyl benzoate (0.5 μ g/mL)