


Analysis of Fluocinolone acetonide

Data No. LA182-0000

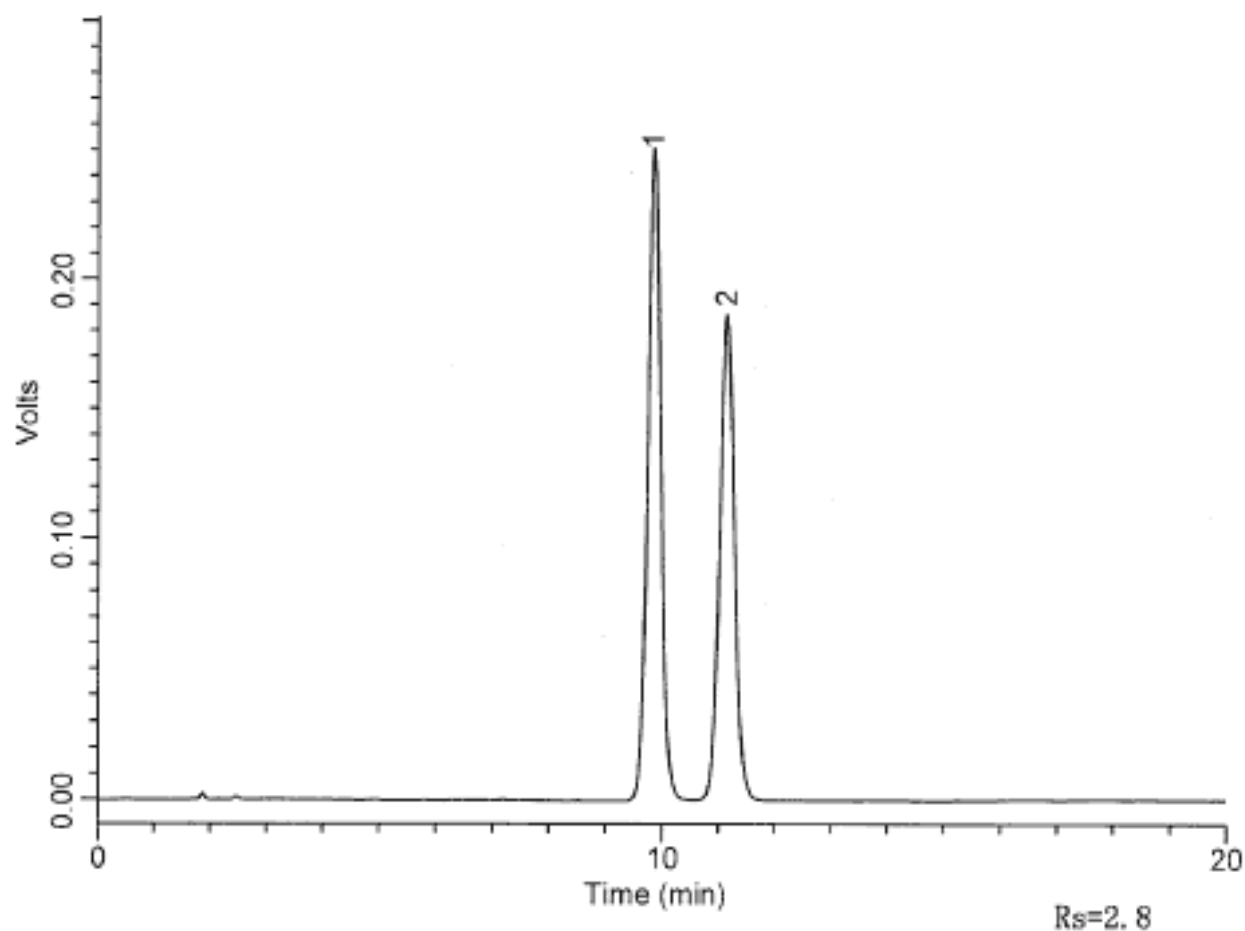

Column : Inertsil SIL-100A
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. Triamcinolone acetonide(1.5mg/mL)
2. Fluocinolone acetonide(1.5mg/mL)

Column Cat. No. :

Eluent : CHCl₃(aq.saturated)/CH₃OH/CH₃COOH
=1000/15/10

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 254 nm

Injection Vol. : 20 μ L

Sample : 1. Triamcinolone acetonide(1.5mg/mL)
2. Fluocinolone acetonide(1.5mg/mL)