


Analysis of Acetaminophen

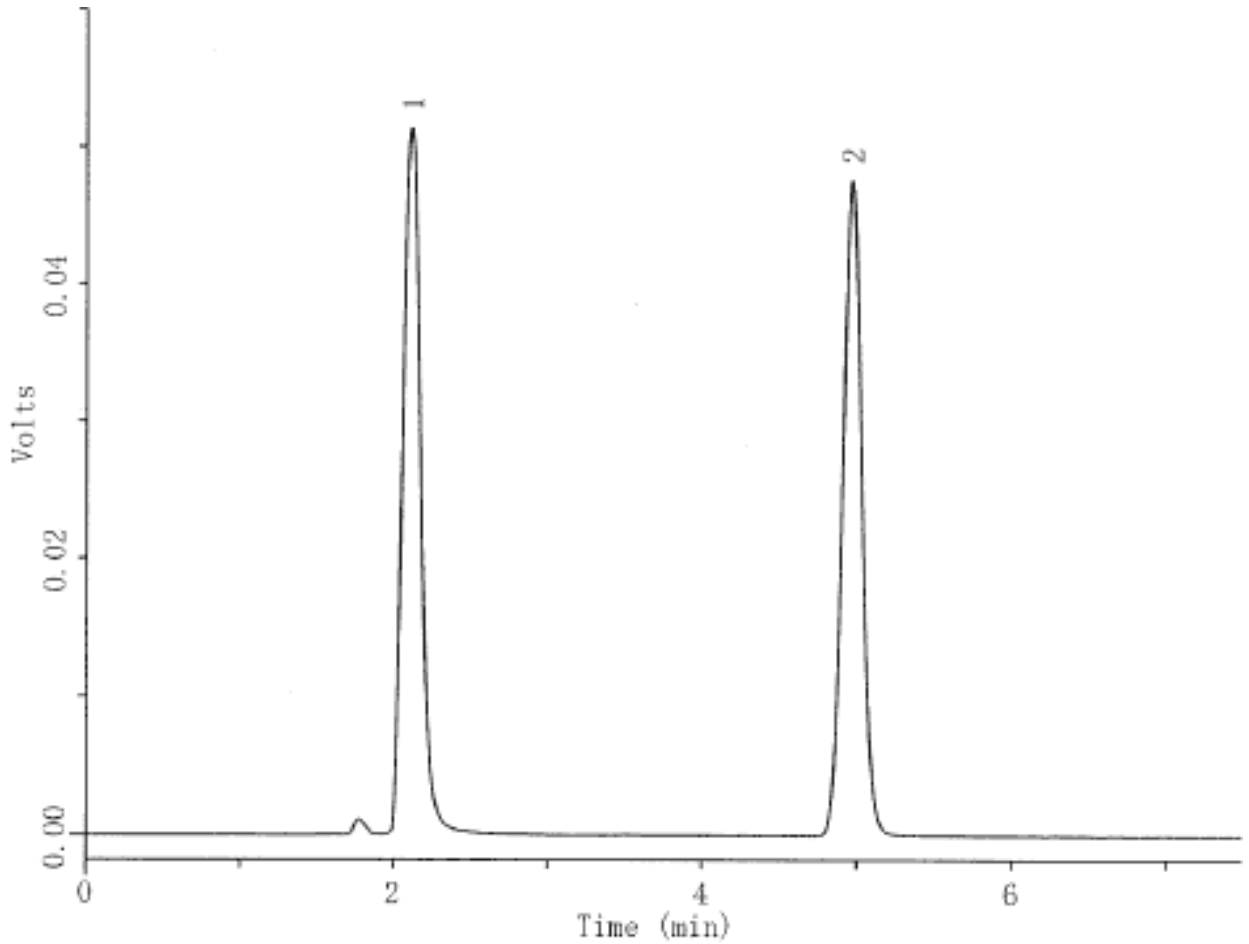

Rs=12.3

Column	:Inertsil ODS-3V 5µm 150 × 4.6mmI.D.	Solutes	: 1. p-Aminophenol(0.02mg/mL) 2. Acetaminophen(0.02mg/mL)
Column Cat. No. :			
Eluent	:CH ₃ CN/50mM KH ₂ PO ₄ =20/80		
Flow rate	:0.9 mL/min.		
Col. Temp.	:40°C		
Detection	:UV 225 nm		
Injection Vol.	:10µL		
Sample	: 1. p-Aminophenol(0.02mg/mL) 2. Acetaminophen(0.02mg/mL)		