


Analysis of Acetaminophen

Data No. LA183-0000

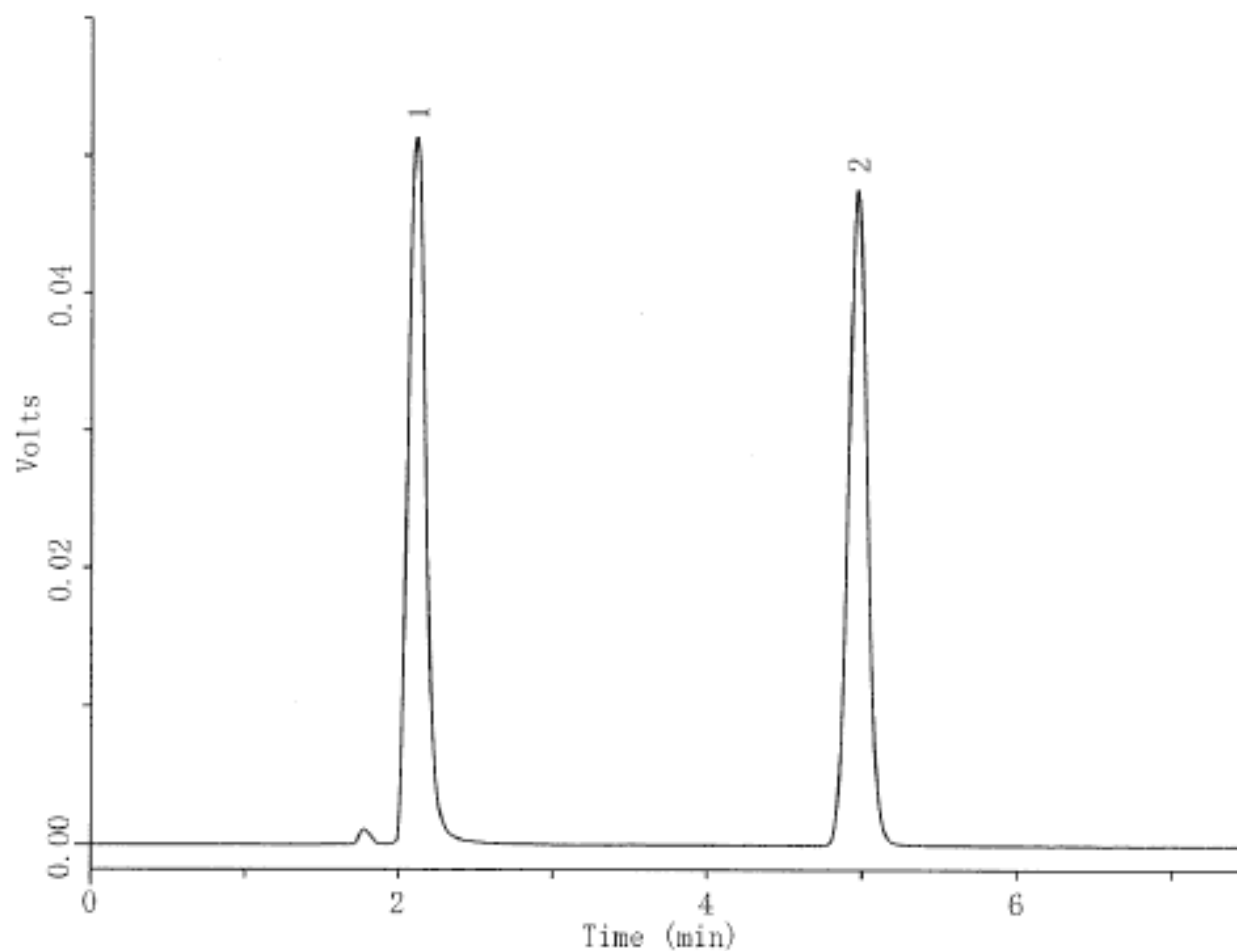

Rs=12.3

Column : Inertsil ODS-3V
5 μ m
150 \times 4.6mm I.D. Solutes : 1. p-Aminophenol(0.02mg/mL)
2. Acetaminophen(0.02mg/mL)

Column Cat. No. :

Eluent : CH₃CN/50mM KH₂PO₄=20/80

Flow rate : 0.9 mL/min.

Col. Temp. : 40°C

Detection : UV 225 nm

Injection Vol. : 10 μ L

Sample : 1. p-Aminophenol(0.02mg/mL)

2. Acetaminophen(0.02mg/mL)