


InertSearch™ for LC

Inertsil® Applications

Analysis of Idoxuridine

Data No. La189-0000

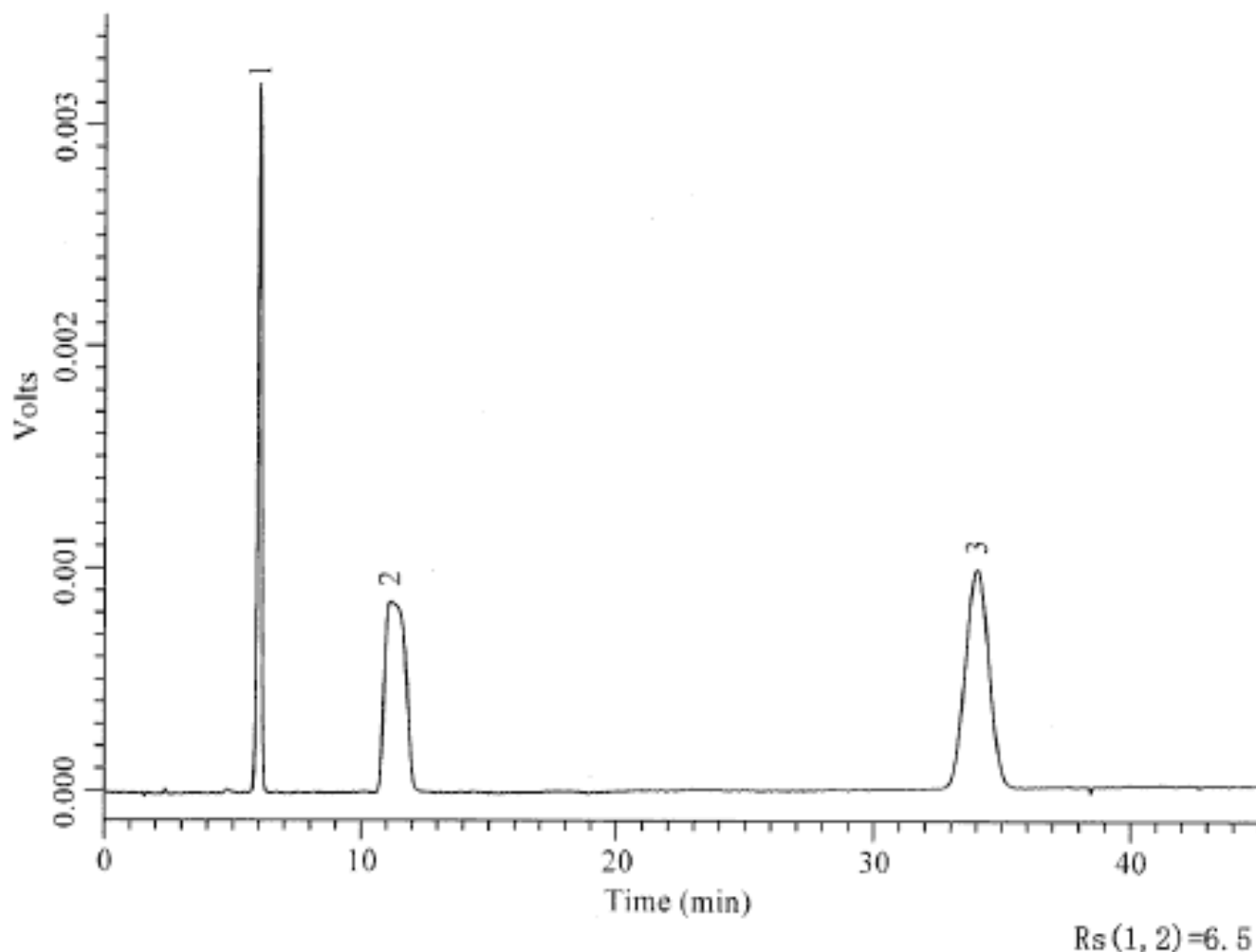

Column : Inertsil ODS-3V
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. 2'-Deoxyuridine (4ng/mL)
2. Iodouracil(12ng/mL)
3. Idoxuridine (20ng/mL)

Column Cat. No. :

Eluent : CH₃OH/H₂O=4/96

Flow rate : 1.3 mL/min.

Col. Temp. : 30°C

Detection : UV 254 nm

Injection Vol. : 10 μ L

Sample : 1. 2'-Deoxyuridine (4ng/mL)
2. Iodouracil(12ng/mL)
3. Idoxuridine (20ng/mL)