


InertSearch™ for LC

Inertsil® Applications

Analysis of Idoxuridine

Data No. LA190-0000

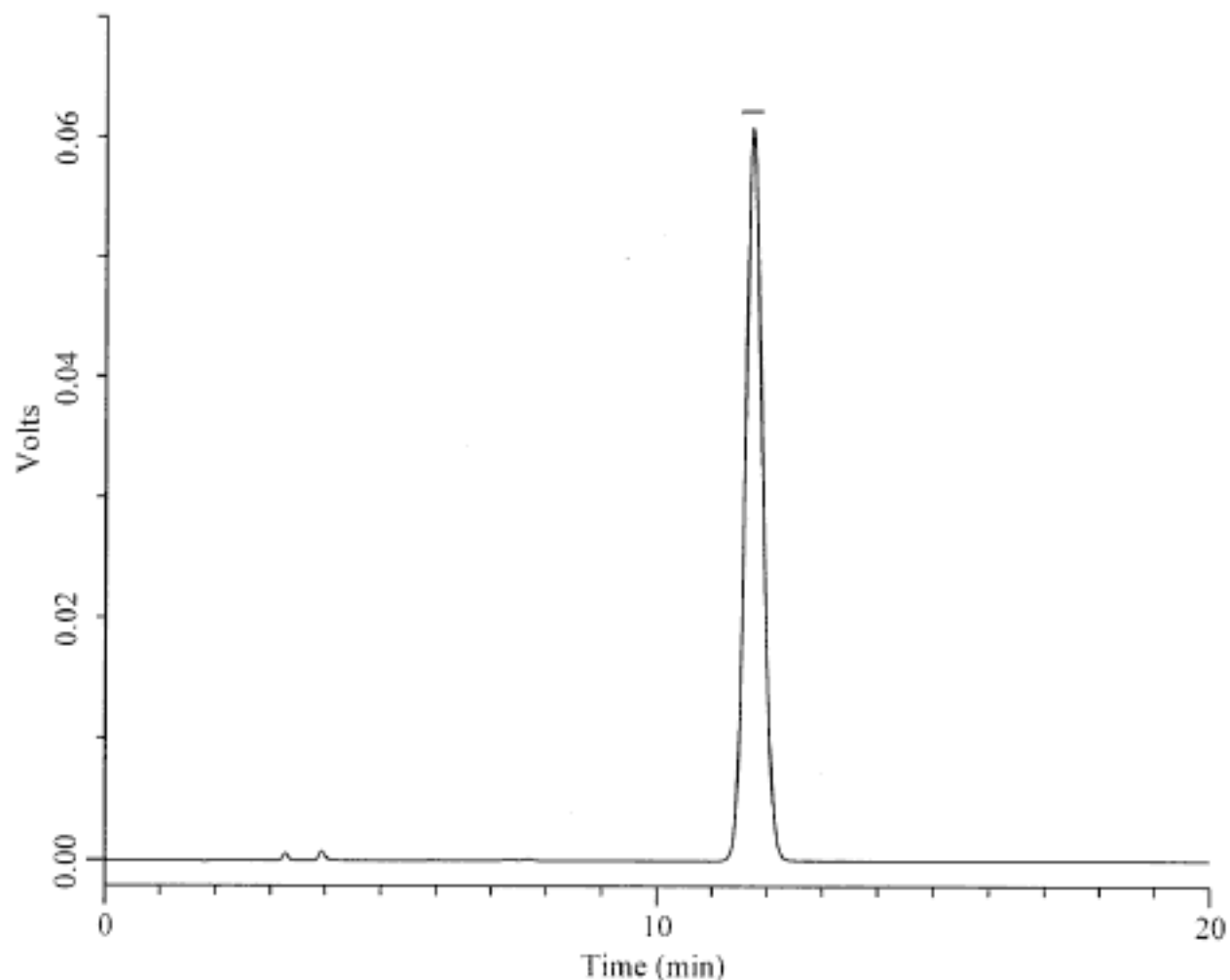

Column : Inertsil ODS-3V
5 μ m
250 \times 4.6mm I.D.

Column Cat. No. :

Eluent : CH₃OH/H₂O=13/87

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 254 nm

Injection Vol. : 10 μ L

Sample : 1. Idoxuridine (0.32mg/mL)

Solutes : 1. Idoxuridine (0.32mg/mL)