


Analysis of Clofibrate

Data No. LA192-0000

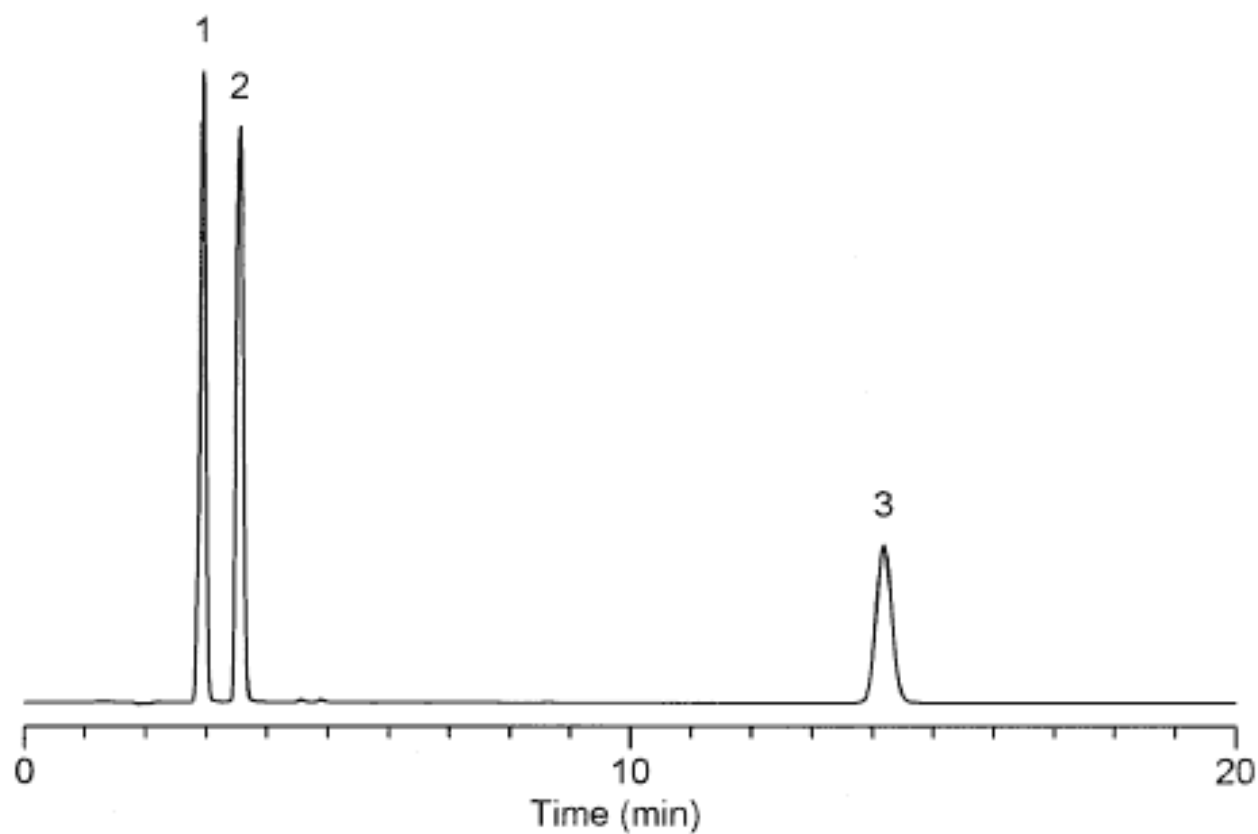

Column :Inertsil ODS-3V
5µm
250 × 4.6mmI.D.

Solutes : 1. p-Ethoxyphenol 2.953
2. p-Chlorophenol 3.563
3. Clofibrate 14.180

Column Cat. No. :

Eluent :acetonitrile / 0.1%phosphoric acid=60/40

Flow rate :1.5 mL/min.

Col. Temp. :25°C

Detection :UV 275 nm

Injection Vol. :10µL

Sample :Clofibrate (1mg/mL)

p-Chlorophenol (0.6mg/mL)

p-Ethoxyphenol (0.6mg/mL)