


Components of Bearberry Leaf

Data No. LA200-0000

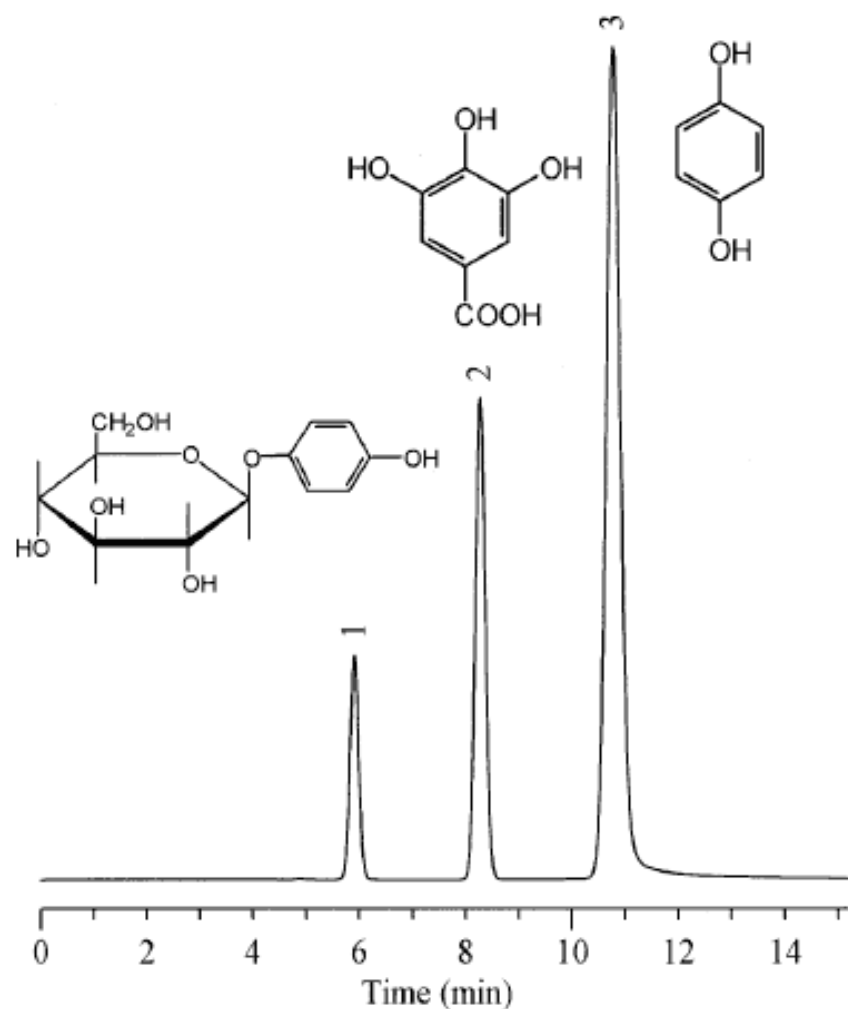

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. Arbutin (4.5 μ g)
2. Gallic acid (4.5 μ g)
3. Hydroquinone (4.5 μ g)

Column Cat. No. :

Eluent : H₂O / MeOH / 0.1N HCl (94 : 5 : 1)

Flow rate : 1.0 mL/min.

Col. Temp. : 25°C

Detection : UV 280 nm

Injection Vol. :

Sample : 1. Arbutin (4.5 μ g)
2. Gallic acid (4.5 μ g)
3. Hydroquinone (4.5 μ g)