
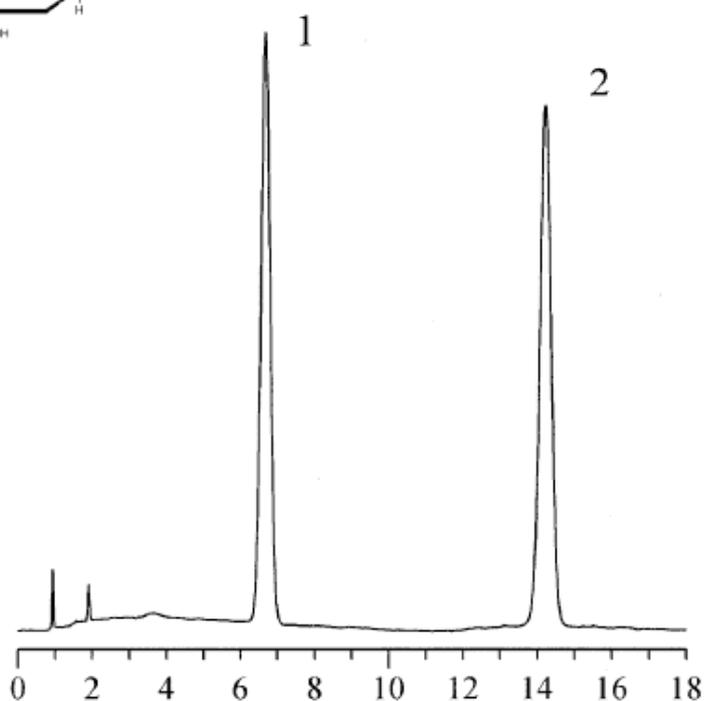


InertSearch™ for LC

Inertsil® Applications

Digitoxin

Data No. LA206-0000

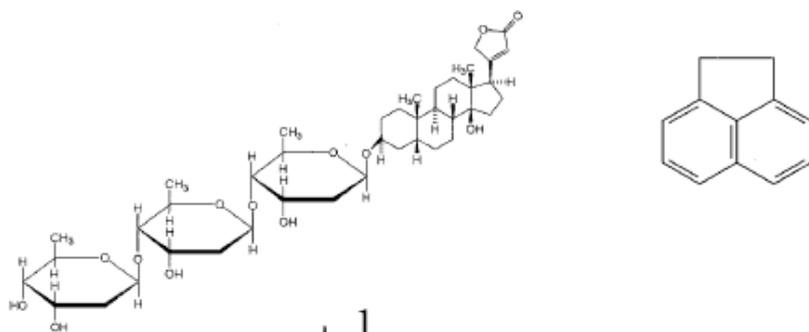

Column : Inertsil ODS-3
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. Digitoxin (0.5 μ g)
2. Accenaphthene (30ng)

Column Cat. No. :

Eluent : MeOH / H₂O (75 : 25) MeOH / H₂O (75 : 25)

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 230 nm

Injection Vol. :

Sample : 1. Digitoxin (0.5 μ g)

2. Accenaphthene (30ng)