


InertSearch™ for LC

Inertsil® Applications

Thiamine

Data No. LA208-0000

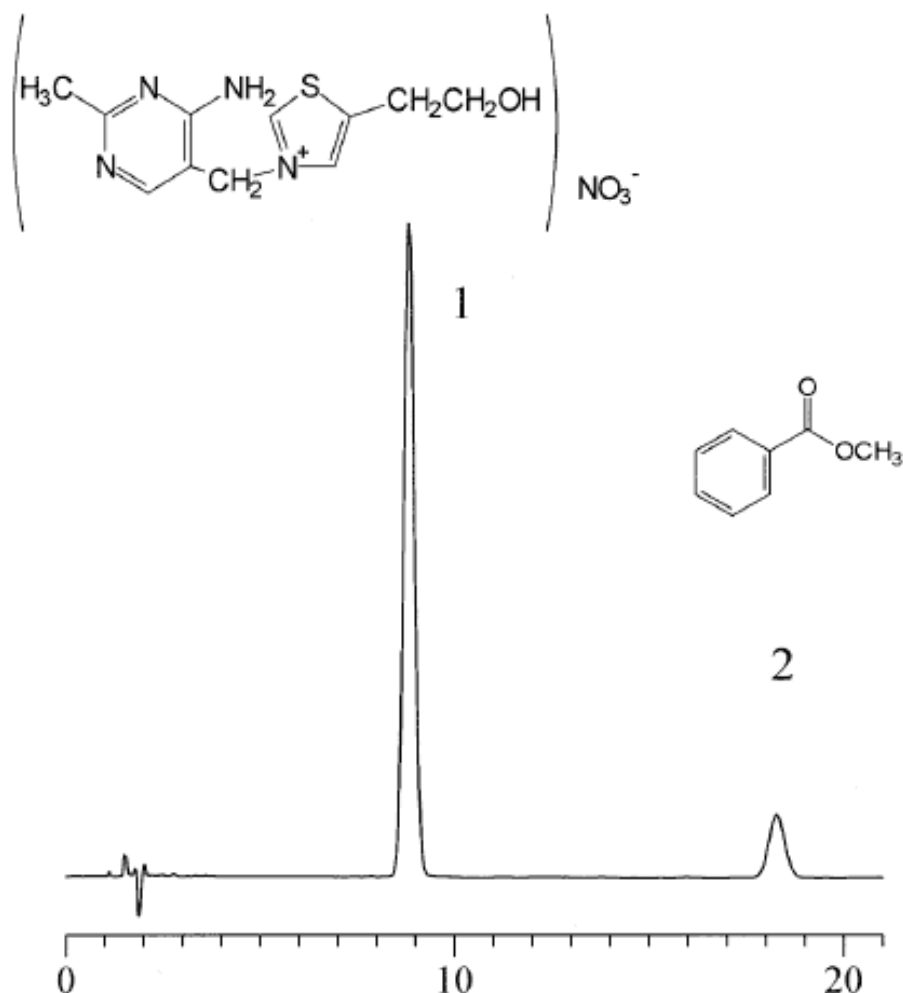

Column : Inertsil ODS-3
5µm
150 × 4.6mm I.D.

Solutes : 1. Thiamine nitrate (2.0µg)
2. Methyl benzoate (0.4 nL)

Column Cat. No. :

Eluent : MeOH / MeCN / 5mM Sodium 1-octanesulfonate in 1%
Acetic acid sol.
= 600 / 240 / 160

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : UV 254 nm

Injection Vol. :

Sample : 1. Thiamine nitrate (2.0µg)
2. Methyl benzoate (0.4 nL)