


Analysis of Cortisone acetate

Data No. LA212-0000

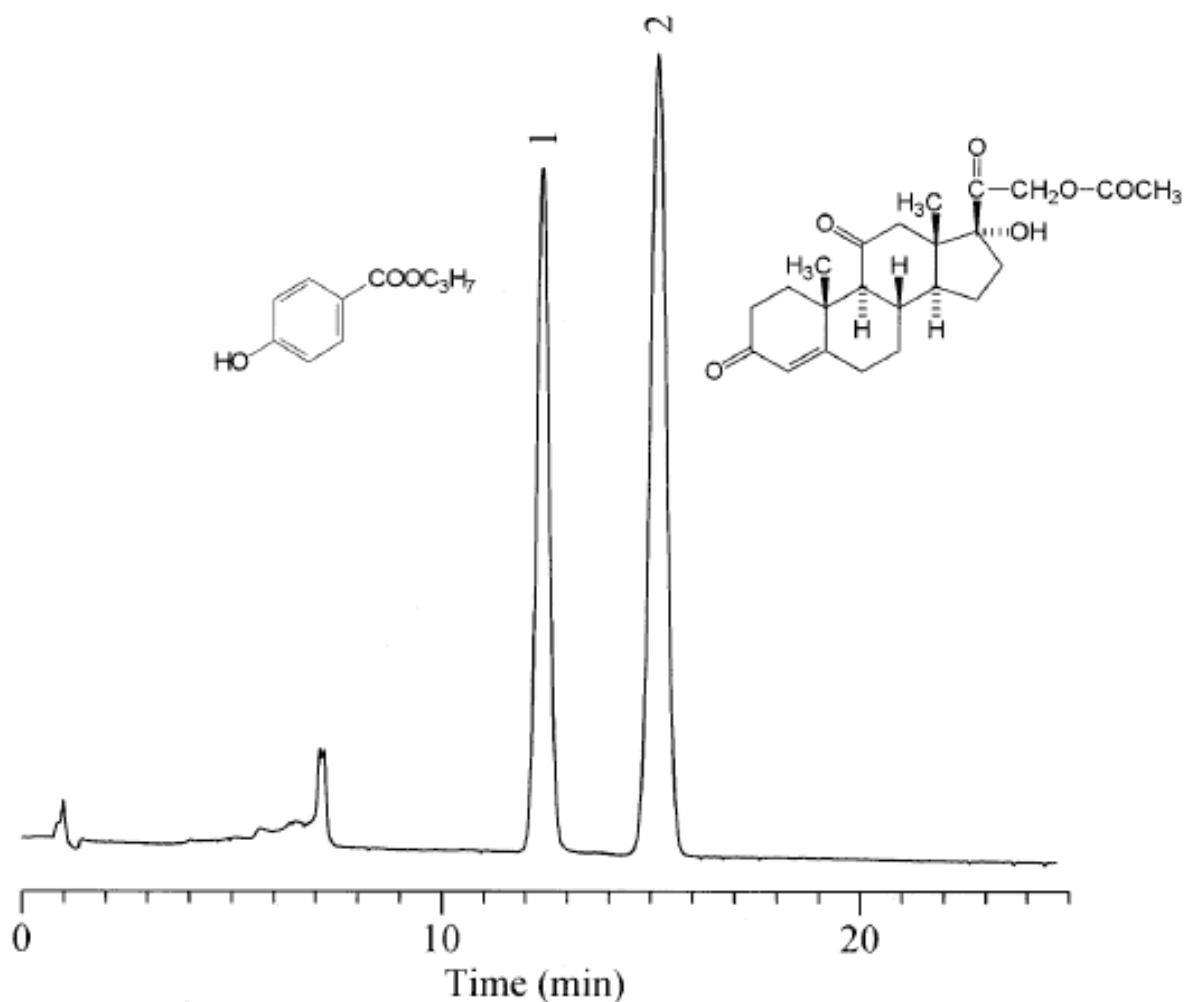

Column : Inertsil ODS-3
5µm
150 × 4.6mm I.D.

Solutes : 1. p-Hydroxy benzoic acid propyl ester 90ng
2. Cortisone acetate 300ng

Column Cat. No. :

Eluent : MeCN / H₂O (35 : 65)

Flow rate : 1.5 mL/min.

Col. Temp. : 25°C

Detection : UV 254 nm

Injection Vol. :

Sample : 1. p-Hydroxybenzoic acid propyl ester 90ng

2. Cortisone acetate 300ng