


Tricyclic antidepressant

Data No. LA213-0000

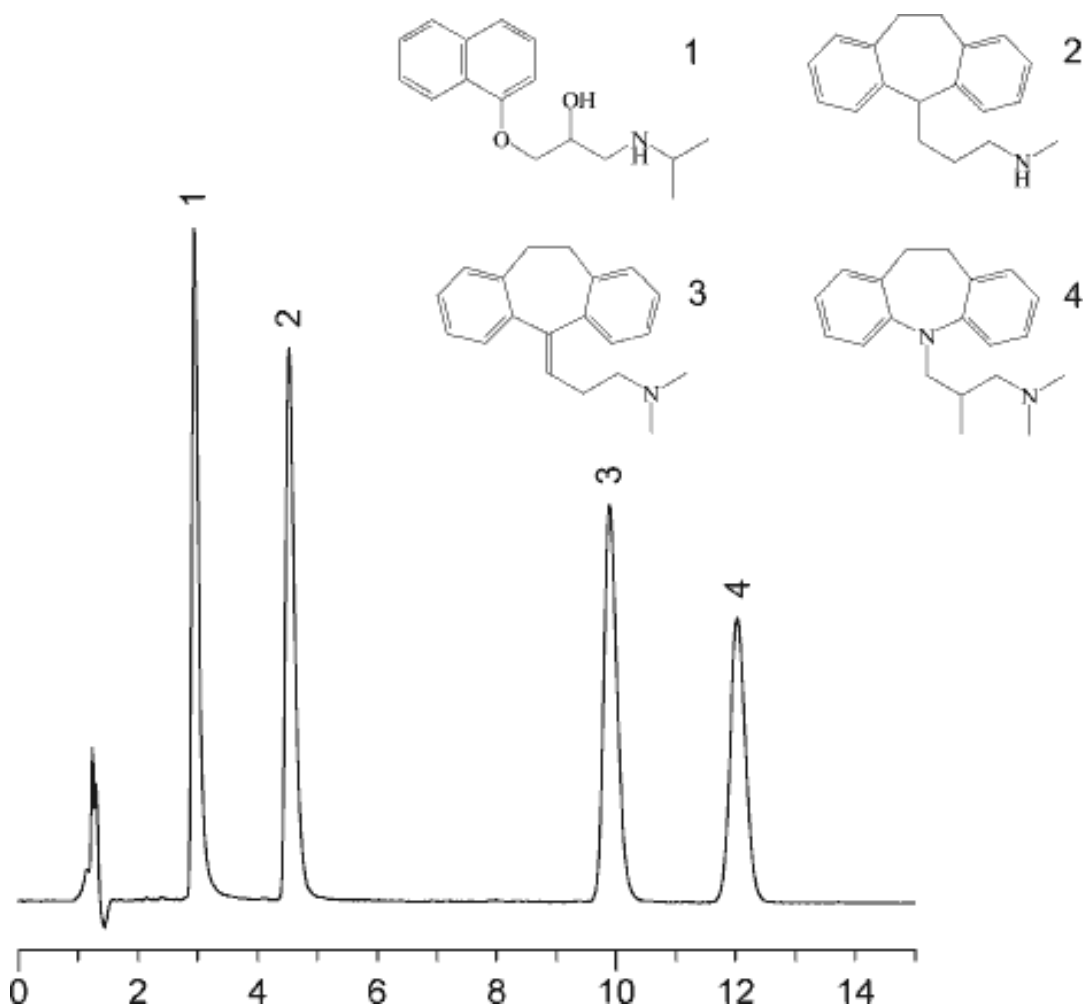

Column :Inertsil ODS-3
5µm
150 × 4.6mmI.D.

Column Cat. No. :

Eluent :20mM K₂HPO₄ (pH7.0) / Acetonitrile (30:70)

Flow rate :1.0 mL/min.

Col. Temp. :40°C

Detection :UV 254 nm

Injection Vol. :3µL

Sample :1. Propranolol (56µg/mL)
2. Nortriptyline (11µg/mL)
3. Amitriptyline (11µg/mL)
4. Trimipramine (11µg/mL)

Solutes : 1. Propranolol (56µg/mL)
2. Nortriptyline (11µg/mL)
3. Amitriptyline (11µg/mL)
4. Trimipramine (11µg/mL)