


Bases of Nucleic acid and Nucleosides

Data No. LA214-0000

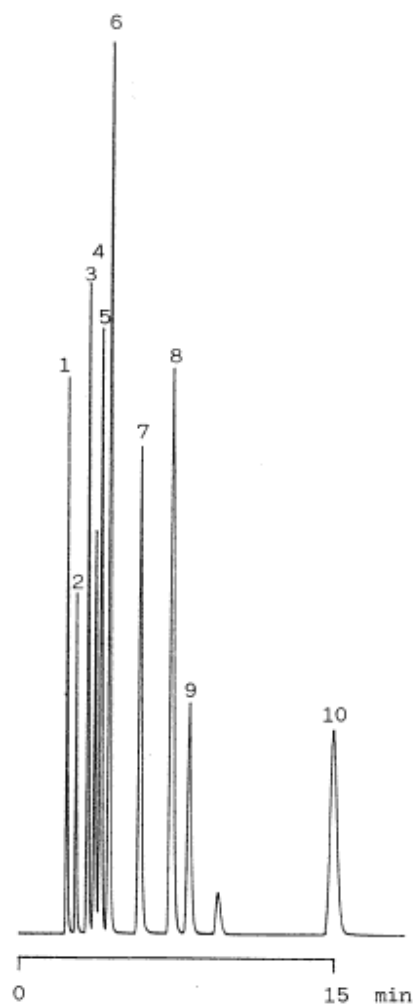

Column : Inertsil ODS-2
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-01101

Eluent : 0.1M H₃PO₄+0.2M NaClO₄(pH2)

Flow rate : 1.0 mL/min.

Col. Temp. : 50°C

Detection : UV 260 nm

Injection Vol. :

Sample :

Solutes : 1. Cytosine
2. Uracil
3. Guanine
4. Adenine
5. Cytidine
6. Uridine
7. Thymine
8. Adenosine
9. Guanosine
10. Thymidine