


InertSearch™ for LC

Inertsil® Applications

Nucleotides

Data No. LA226-0000

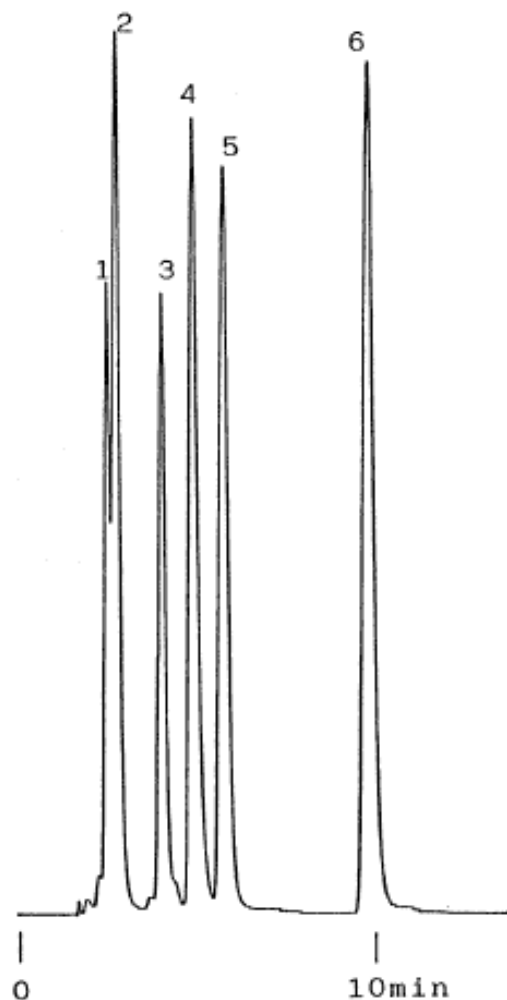

Column :Inertsil ODS-2
5 μ m
150 \times 0.7 mmI.D.

Column Cat. No. :5020-13073

Eluent :0.1M Na₂HPO₄(pH6 by H₃PO₄)

Flow rate :30 μ L/min.

Col. Temp. :

Detection :UV 260 nm

Injection Vol. :

Sample :

Solutes : 1. GTP
2. GDP
3. GMP
4. ATP
5. ADP
6. AMP