


Analysis of Adrenocortical hormones

Data No. LA241-0000

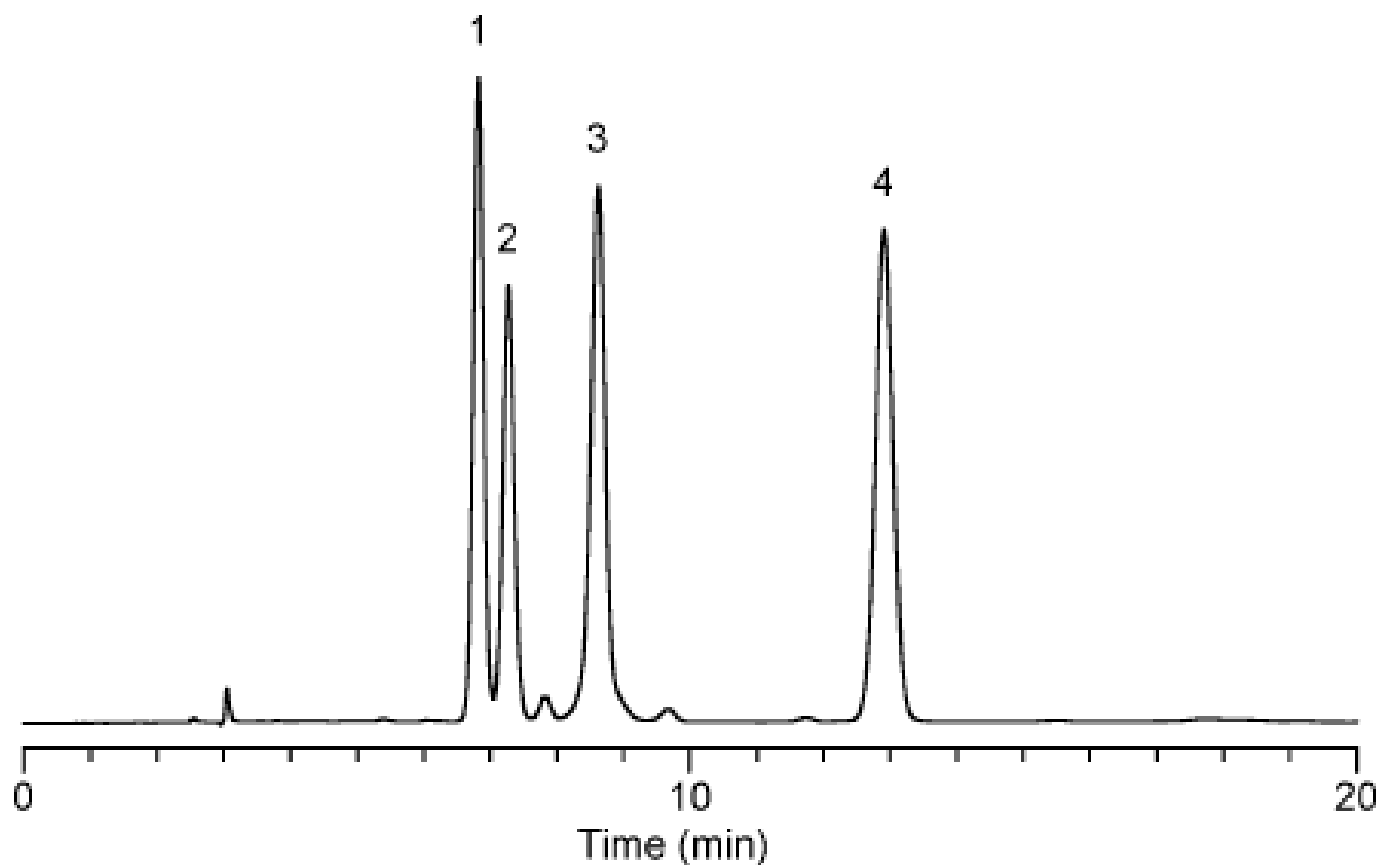

Column	: Inertsil ODS-3V 5µm 150 × 4.6mmI.D.	Solutes	: 1. Prednisone 2. Cortisone 3. Prednisolone 4. Corticosterone
Column Cat. No.	:		
Eluent	: Methanol/Water=60/40		
Flow rate	: 1.0mL/min		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	: 3µL		
Sample	: Prednisone (0.1mg/mL) Cortisone (0.1mg/mL) Prednisolone (0.1mg/mL) Corticosterone (0.1mg/mL)		