


Analysis of Gamma-Oryzanol

Data No. LA248-0000

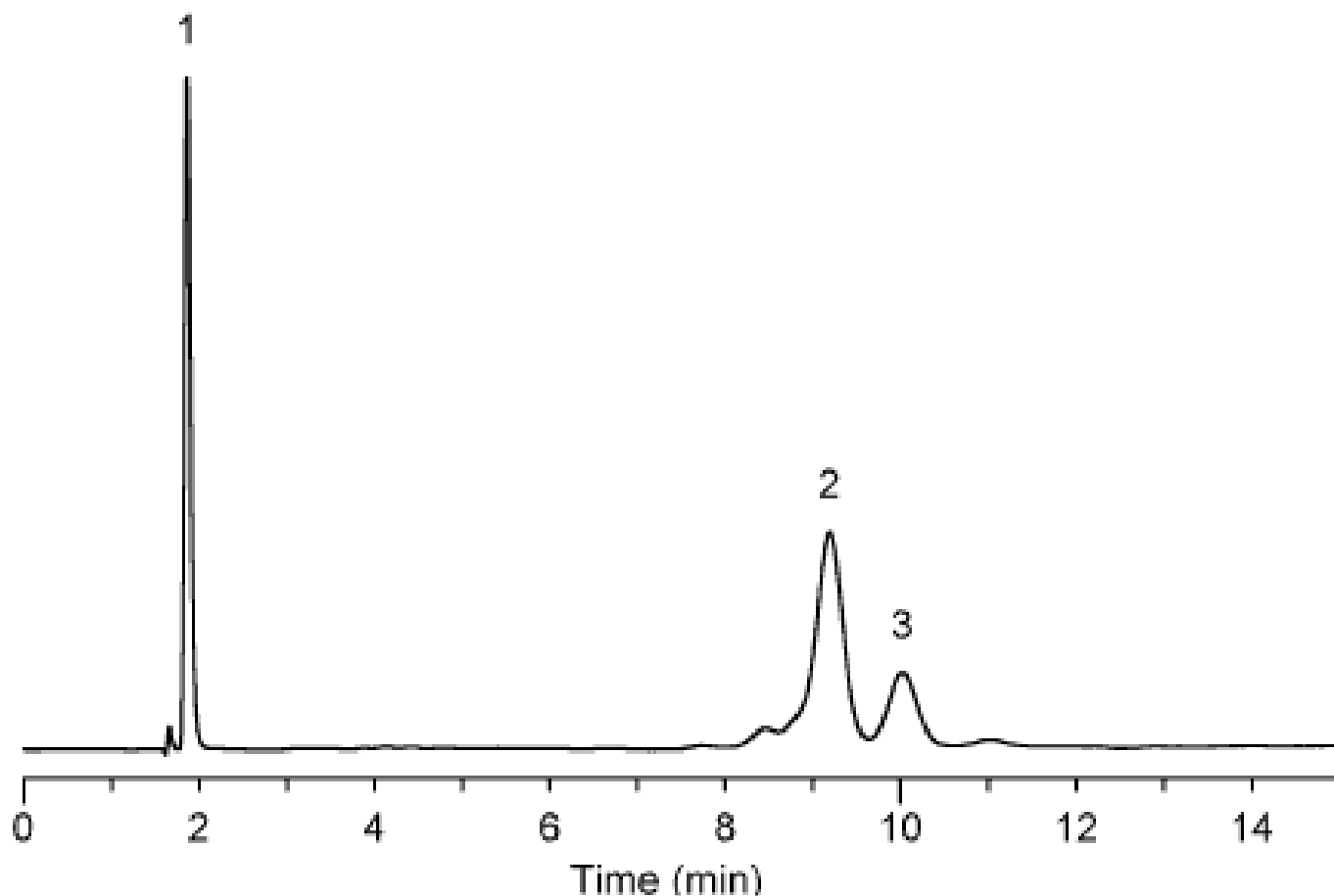

Column : Inertsil ODS-3V
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. vanillin
2, γ -oryzanol-1
3. γ -oryzanol-2

Column Cat. No. :

Eluent : Methanol:2-Propanol:Water=70:30:1

Flow rate : 1.0mL/min

Col. Temp. : 40°C

Detection : UV 320 nm

Injection Vol. : 1 μ L

Sample : γ -Oryzanol