


InertSearch™ for LC

Inertsil® Applications

Analysis of Gamma-Oryzanol

Data No. LA249-0000

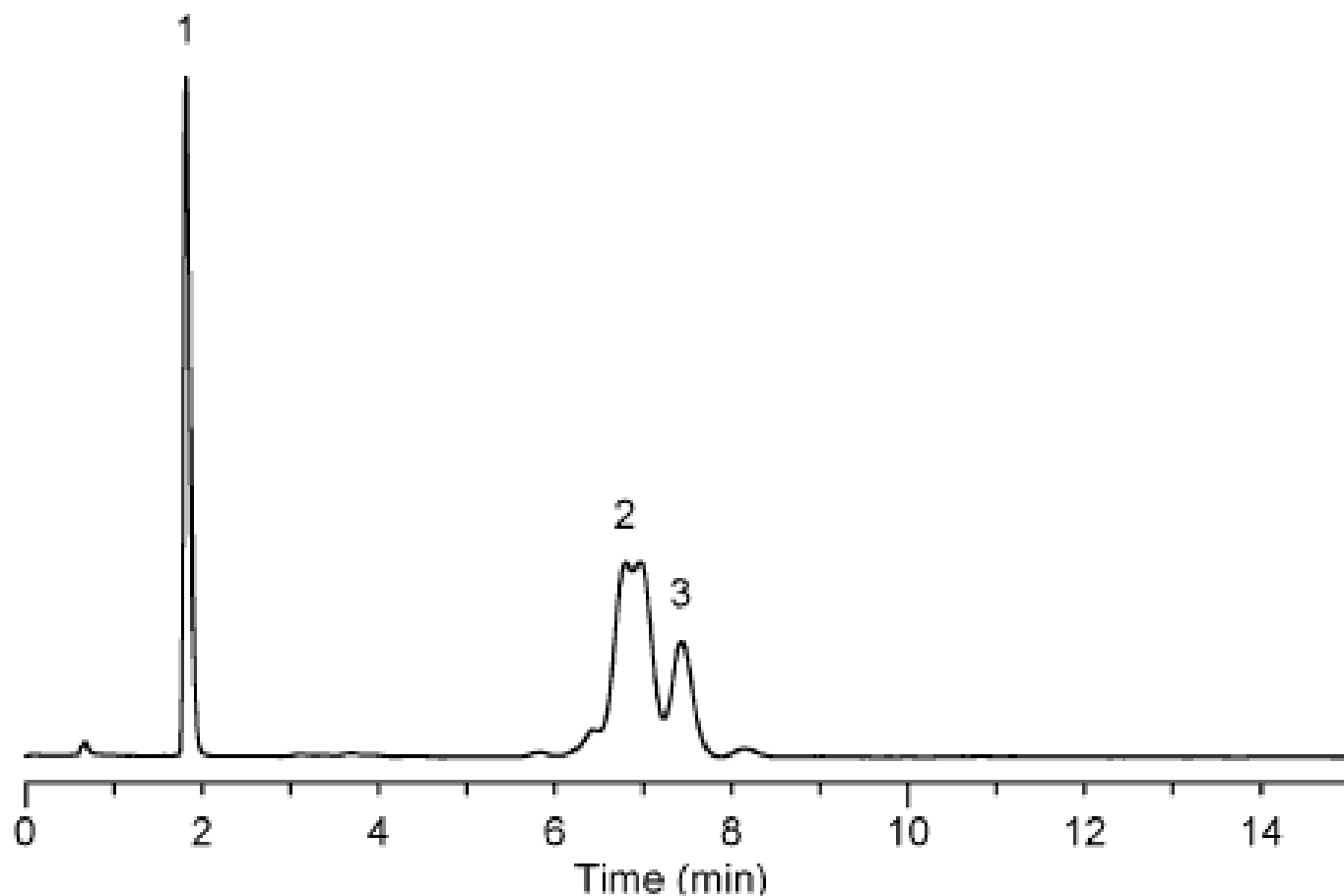

Column : Inertsil C8-3
5 μ m
150 \times 4.6mm I.D.

Solutes : 1. vanillin (50ug/mL)
2. γ -oryzanol (mixture 100ug/mL)

Column Cat. No. : 5020-01900

Eluent : Methanol: Water=95:5

Flow rate : 1.0mL/min

Col. Temp. : 40°C

Detection : UV 320 nm

Injection Vol. : 1 μ L

Sample : γ -Oryzanol