


## Analysis of Steroid Hormones

Data No. LA251-0000

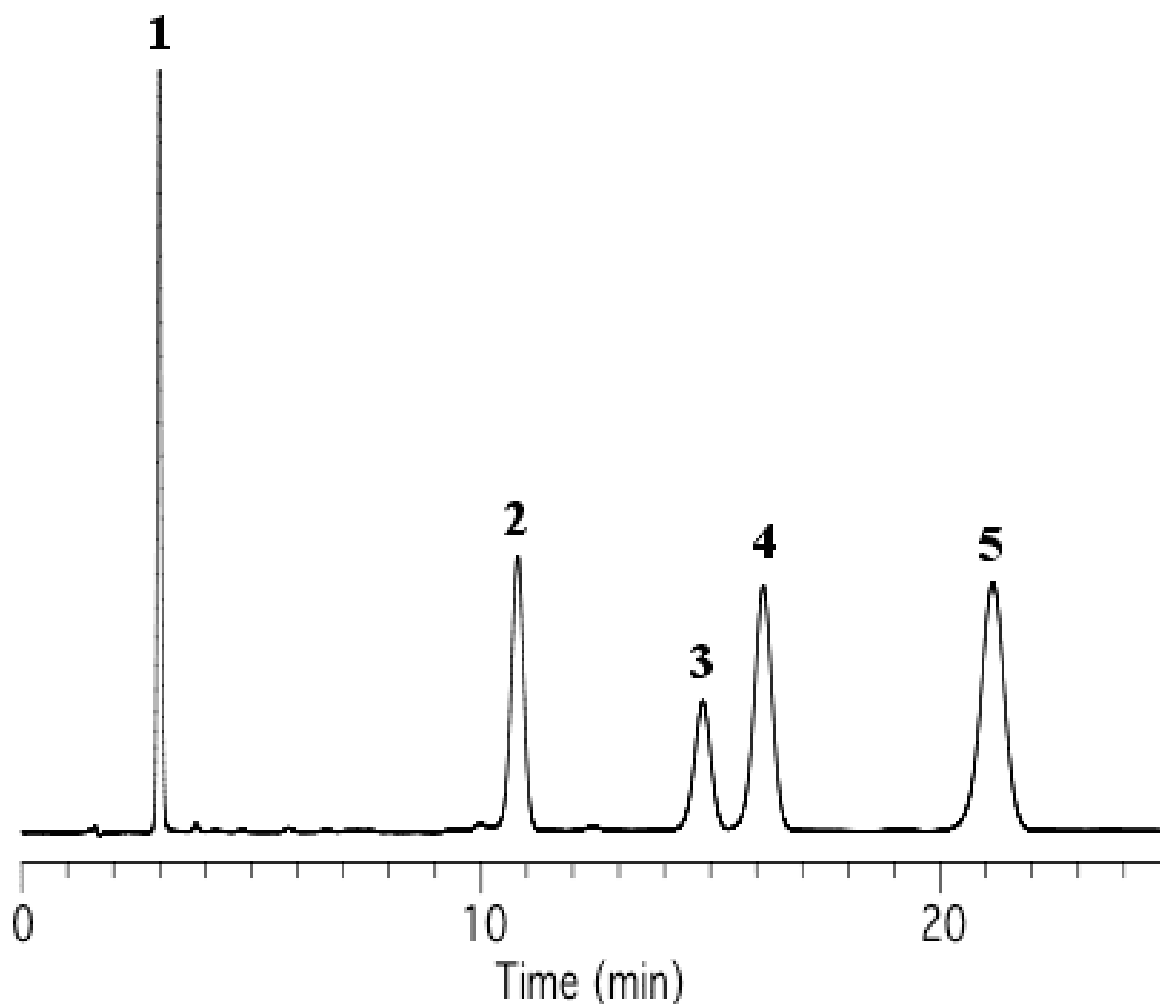

Column	: Inertsil Diol 5 µm 150 x 4.6 mm I.D.	Solutes	: 1. Progesterone 2. Corticosterone 3. Cortisone 4. Prednisone 5. Prednisolone
Column Cat. No.	: 5020-05645		
Eluent	: A) n-Hexane B) Ethanol A/B = 85/15, w/w		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		