


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA256-0000

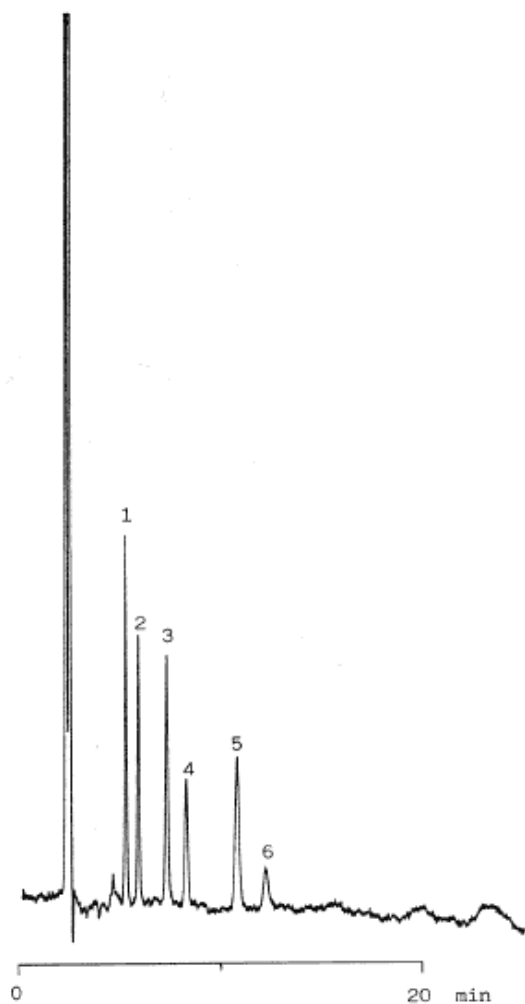

Column :Inertsil SIL-150A
5µm

Column Cat. No. :5020-01101

Eluent :n-Hexane/C₂H₅OH=99/1

Flow rate :1.0 mL/min.

Col. Temp. :40°C

Detection :UV 230nm

Injection Vol. :

Sample :

Solutes : 1. Prometryne
2. Terbutryne
3. Ametryne
4. Atrazine
5. Simetryne
6. Sinazine