


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA262-0000

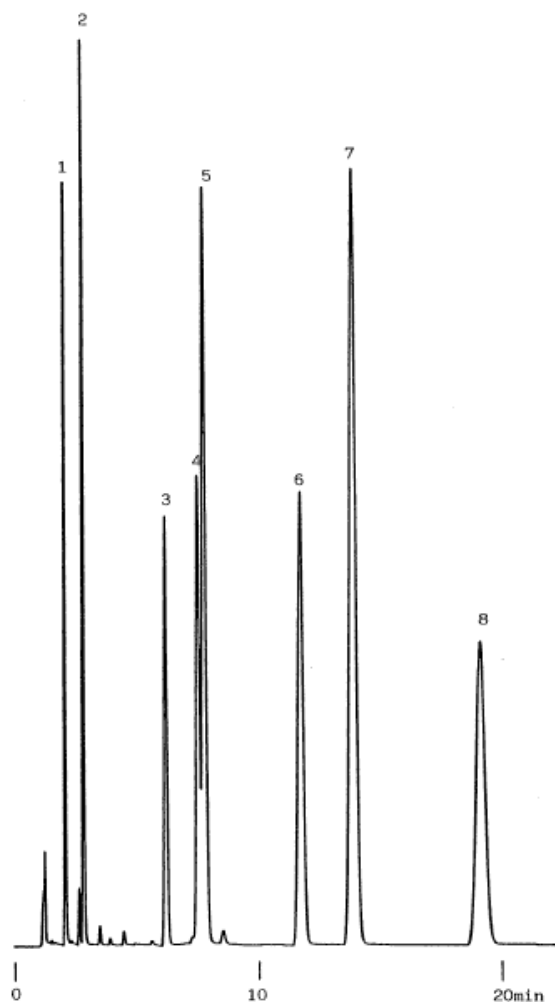

Column : Inertsil ODS-2
5 μ m
150 \times 4.6mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃CN/H₂O=55/45

Flow rate : 1.0 mL/min.

Col. Temp. : 50°C

Detection : UV 235 nm

Injection Vol. :

Sample :

Solutes : 1. Dimethoate
2. DDVP
3. Methyl parathion
4. Malathion
5. Sumithion
6. Parathion
7. Diazinon
8. EPN
each 0.1%