


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA266-0000

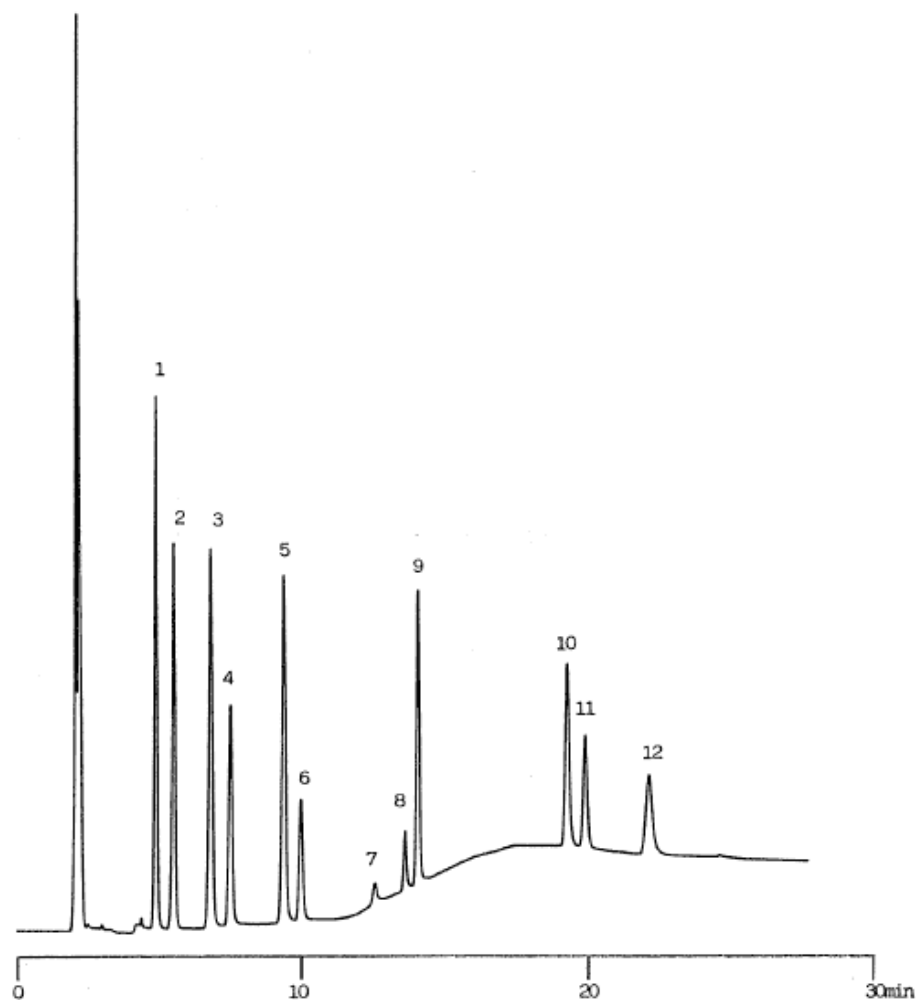

Column	:Inertsil SIL-150A 5µm	Solutes	: 1. Prometryne 2. Terbutryne 3. Ametryne 4. Atrazine 5. Simetryne 6. Sinazine 7. Propoxur 8. Carbofuran 9. Carbaryl 10. Fenuron 11. Monuron 12. Methomyl 13. each 40ng
Column Cat. No.	:5020-01011		
Eluent	:A:n-Hexane B:C ₂ H ₅ OH A/B 99/1→(6min)→98/2 →(5min)?90/10		
Flow rate	:1.0 mL/min.		
Col. Temp.	:40°C		
Detection	:UV 230 nm		
Injection Vol.	:		
Sample	:		