


InertSearch™ for LC

Inertsil® Applications

Pesticides

Data No. LA270-0000

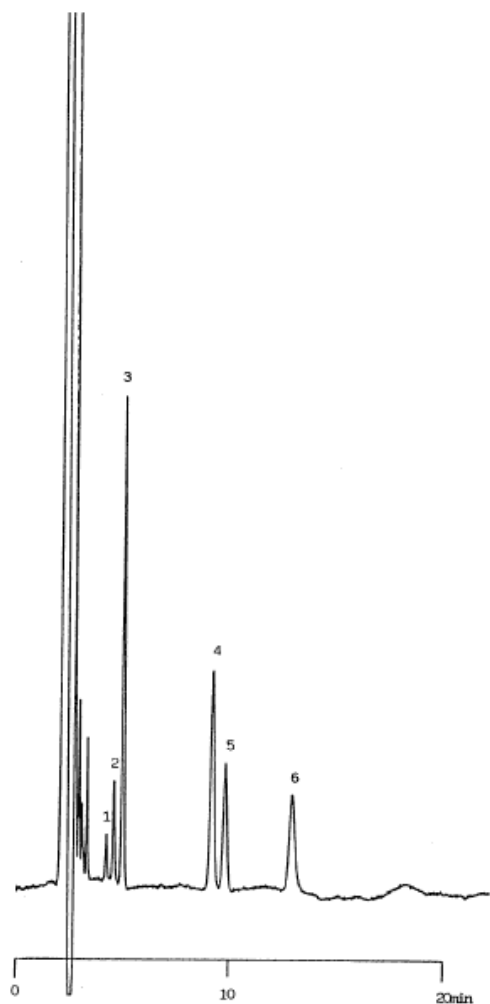

Column :Inertsil SIL-150A
5µm

Column Cat. No. :5020-01011

Eluent :n-Hexane/C₂H₅OH=96/4

Flow rate :1.0 mL/min.

Col. Temp. :40°C

Detection :UV 230 nm

Injection Vol. :

Sample :

Solutes : 1. Propoxur
2. Carbofuran
3. Carbaryl
4. Fenuron
5. Monuron
6. Methomyl